



ease@T

SOFT, SAFE, SYNCED BREAST HEALTH


SCAN CENTRE PERFORMING PROVIDER USER MANUAL

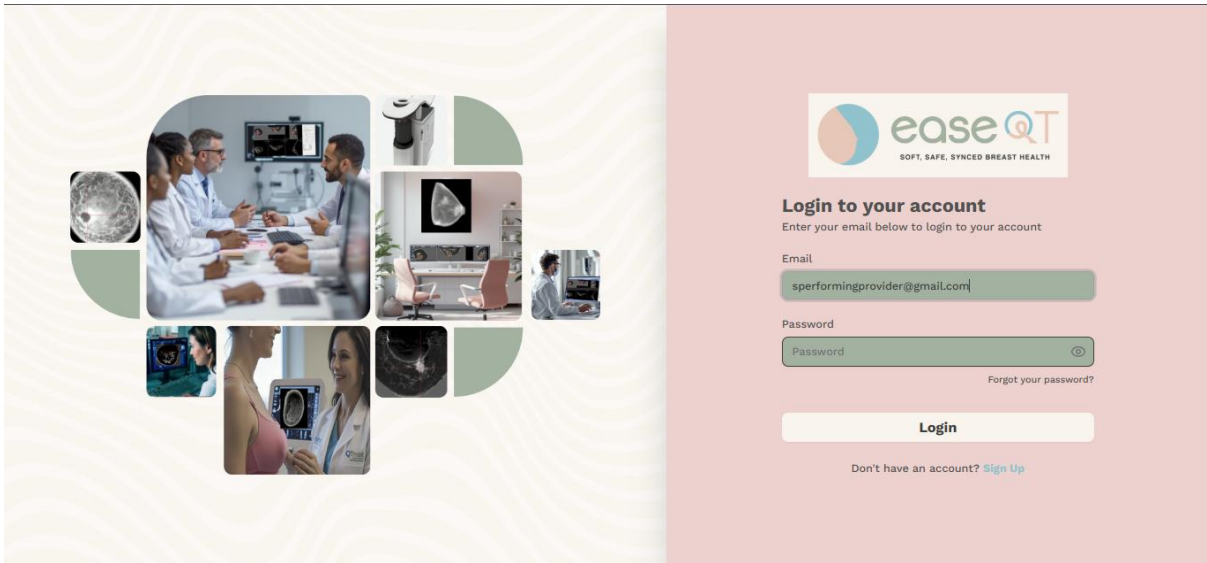

The performing provider can sign in using email and password provided by their concern scan centre.

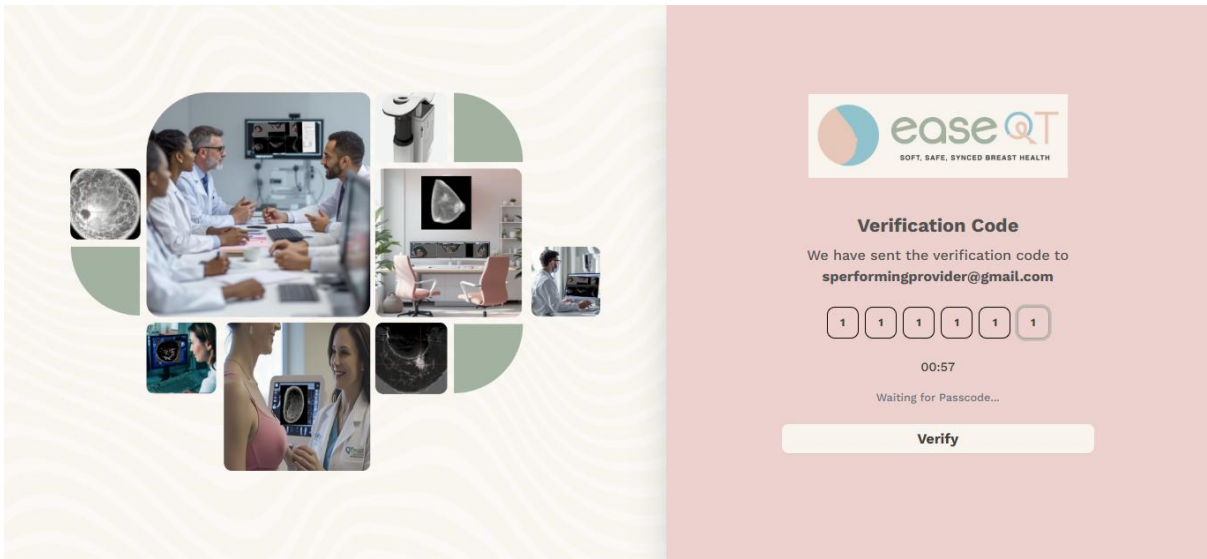


A verification code will be generated and sent for verification.

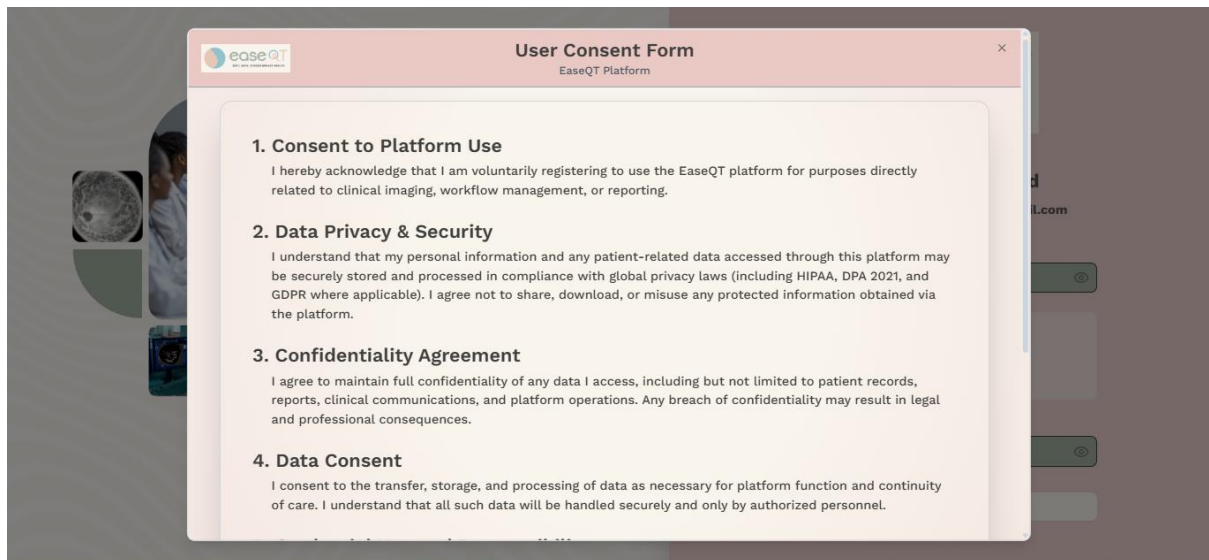

The screenshot shows a web browser window with the title "User Consent Form" and the subtitle "EaseQT Platform". The form is divided into four sections:

- 1. Consent to Platform Use**
I hereby acknowledge that I am voluntarily registering to use the EaseQT platform for purposes directly related to clinical imaging, workflow management, or reporting.
- 2. Data Privacy & Security**
I understand that my personal information and any patient-related data accessed through this platform may be securely stored and processed in compliance with global privacy laws (including HIPAA, DPA 2021, and GDPR where applicable). I agree not to share, download, or misuse any protected information obtained via the platform.
- 3. Confidentiality Agreement**
I agree to maintain full confidentiality of any data I access, including but not limited to patient records, reports, clinical communications, and platform operations. Any breach of confidentiality may result in legal and professional consequences.
- 4. Data Consent**
I consent to the transfer, storage, and processing of data as necessary for platform function and continuity of care. I understand that all such data will be handled securely and only by authorized personnel.


The consent form must be filled by the performing provider.

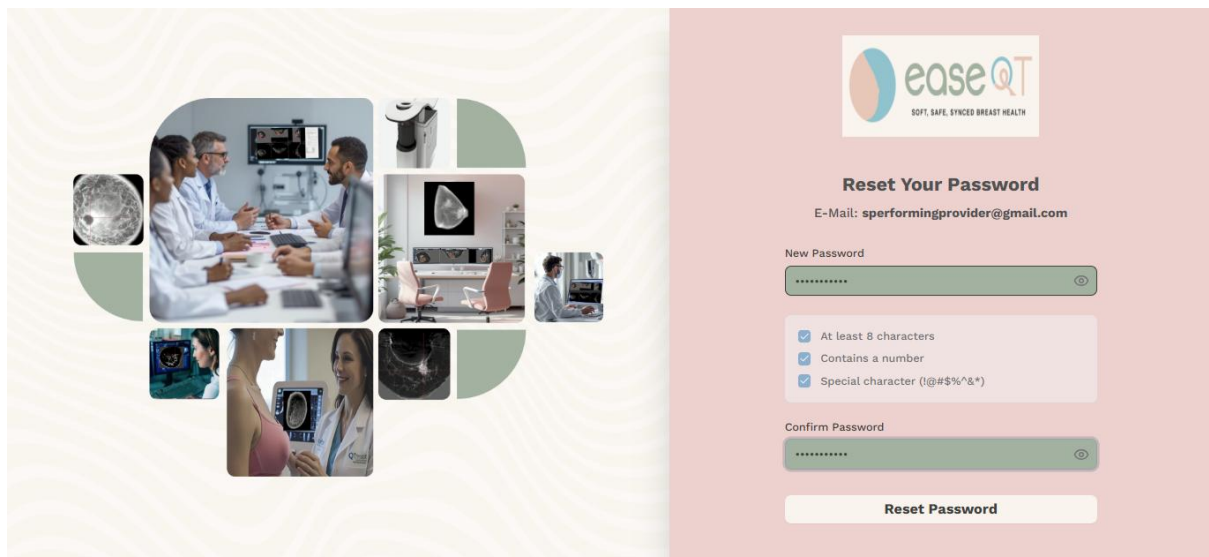

The screenshot shows a web browser window with the title "Reset Your Password" and the subtitle "EaseQT Platform". The form is divided into two main sections:

- Left Section:** A collage of images showing medical professionals in a clinical setting, including a doctor, a nurse, and a patient, along with medical equipment like an ultrasound machine and a monitor.
- Right Section:** A form for resetting the password. It includes the EaseQT logo, the text "Reset Your Password", the email address "E-Mail: `sperformingprovider@gmail.com`", and two password input fields labeled "New Password" and "Confirm Password". Below the input fields are three checkboxes for password requirements: "At least 8 characters", "Contains a number", and "Special character (!@#%&^*)". A "Reset Password" button is at the bottom.

The password must be reset for the security purpose.

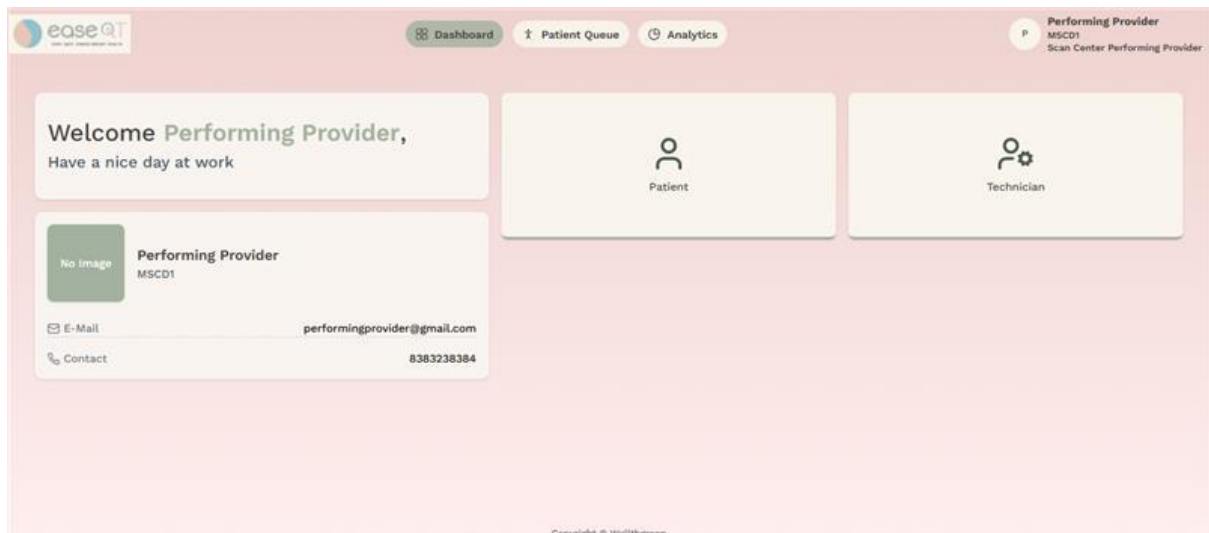

In dashboard, the scan centre performing provider can view patient brochure, patient consent technician guidelines, technician consent.

PATIENT QUEUE

The Patient Queue interface includes a top navigation bar with 'Dashboard', 'Patient Queue', and 'Analytics' tabs. A user profile card on the right displays the provider's name, email (CMSCPP1), and contact number (Scan Center Performing Provider). Below the navigation bar are buttons for 'Download Dicom', 'Download Final Reports', and 'Clear All Filters'. A search bar is also present. The main table lists patient records with columns for Date, Scan Centre, Patient Name, Form, Tech Form, DICOM (L, R), Old Report, QT Report, Report Status, Assign, Pending Remarks, and Remarks. The table is filtered for 'MSC' and 'Testing12'. The bottom of the interface shows pagination controls (Rows per page: 50, Page 1 of 1) and a copyright notice 'Copyright © Wellthgreen'.

	Date	Scan Centre	Patient Name	Form	Tech Form	DICOM L R	Old Report	QT Report	Report Status	Assign	Pending Remarks	Remarks
<input type="checkbox"/>	Sep 17, 2025	MSC	Testing12	Dc - View	View	(1) -		View	SIGNED OFF (A)	<input type="text"/>	No remarks...	+
<input type="checkbox"/>	Sep 18, 2025	MSC	Testing12	S - View	View	(1) -		View	REVIEWED 2	Assign <input type="text"/>	R1 R2 stuck	+
<input type="checkbox"/>	Sep 22, 2025	MSC	Testing12	S - View	View	(1) -		View	SIGNED OFF (A)	Assign <input type="text"/>	No remarks...	+
<input type="checkbox"/>	Sep 18, 2025	MSC	sofia	S - View	View	(1) -		View	DRAFT	Assign <input type="text"/>	No remarks...	+

In the patient queue, the Scan centre performing provider can view the patient list, patient forms, technician view and DICOM files. They can fill out reports, check the report status, assign the case to another user and add pending remarks if needed and add pending remark and can view previous remark also. Click on **View** to proceed with filling out the report. And they can multi select and download DICOM and Final reports. They can upload and download patient's old reports.

Other Old reports

S.No	Report Title	File
1	Thermogram	View

Time stamp

Created by	Start	End	Status
Wellthgreen Performing Provider	19 Sept 2025 21:42	19 Sept 2025 21:57	Saved Changes
Wellthgreen Performing Provider	19 Sept 2025 21:57	19 Sept 2025 22:00	Saved Changes
MSC Performing Provider	19 Sept 2025 22:00	19 Sept 2025 22:01	Saved Changes
MSC Performing Provider	19 Sept 2025 22:01	-	-

Status saver

Load Template

Reset to Default

Predraft

Draft

Reviewed 1 Correct

Reviewed 1 Edit

Reviewed 2 Correct

Reviewed 2 Edit

Insert Signature

Sign Off

Normal

NAME

sofa

DOB

September 2, 1997 (28 y)

GENDER

F

SCAN CENTER

MSC, Bangalore

USERID

P100016

DATE OF VISIT

Sep 18, 2025

QT ULTRASOUND BREAST IMAGING

Indication: Screening : Routine Annual Checkup.

History:

Reason for having this QT scan (Patient Stated): ssss.

Current breast symptoms:

- Lump or thickening in the left breast at 6 o'clock since 2 week, size of a pea, right breast at 2 o'clock since 2 week, size of a pea.
- Skin changes in the left breast at 8 o'clock since 2 week, dimpling is noted, right breast at 3 o'clock since 2 week, dimpling is noted.

TECHNIQUE: Transmission and reflection multiplanar 3-dimensional ultrasound imaging of both breasts was performed using the QT Ultrasound Series 2000 Model-A scanner. Breast density was determined using the Quantitative Breast Density calculator. Images were reviewed in the QTViewer v2.6.2. The nipple-areolar complex, skin, Cooper's ligaments, breast fat distribution, penetrating arteries and veins, glandular and ductal tissues were evaluated. Images were reviewed in coronal, transaxial and sagittal planes.

The scan centre performing provider can review the report and click “Sign Off” to complete it. If the performing provider needs to edit anything in the final report, after completing the edits, they can insert their signature by clicking “Insert Signature. If admin allows the provider can view general, left, right.

Other Old reports

S.No	Report Title	File
1	Thermogram	View

Time stamp

Created by	Start	End	Status
Wellthgreen Performing Provider	19 Sept 2025 21:42	19 Sept 2025 21:57	Saved Changes
Wellthgreen Performing Provider	19 Sept 2025 21:57	19 Sept 2025 22:00	Saved Changes
MSC Performing Provider	19 Sept 2025 22:00	19 Sept 2025 22:01	Saved Changes
MSC Performing Provider	19 Sept 2025 22:01	-	-

Status saver

Load Template

Reset to Default

Predraft

Draft

Reviewed 1 Correct

Reviewed 1 Edit

Reviewed 2 Correct

Reviewed 2 Edit

Insert Signature

Sign Off

Normal

NAME

sofa

DOB

September 2, 1997 (28 y)

GENDER

F

SCAN CENTER

MSC, Bangalore

USERID

P100016

DATE OF VISIT

Sep 18, 2025

QT ULTRASOUND BREAST IMAGING

Indication: Screening : Routine Annual Checkup.

History:

Reason for having this QT scan (Patient Stated): ssss.

Current breast symptoms:

- Lump or thickening in the left breast at 6 o'clock since 2 week, size of a pea, right breast at 2 o'clock since 2 week, size of a pea.
- Skin changes in the left breast at 8 o'clock since 2 week, dimpling is noted, right breast at 3 o'clock since 2 week, dimpling is noted.

TECHNIQUE: Transmission and reflection multiplanar 3-dimensional ultrasound imaging of both breasts was performed using the QT Ultrasound Series 2000 Model-A scanner. Breast density was determined using the Quantitative Breast Density calculator. Images were reviewed in the QTViewer v2.6.2. The nipple-areolar complex, skin, Cooper's ligaments, breast fat distribution, penetrating arteries and veins, glandular and ductal tissues were evaluated. Images were reviewed in coronal, transaxial and sagittal planes.

There is an option called "Reset to Default." If the report was done incorrectly, they can use this option to revert it to the default state.

Patient Form

Technician Form

Dicom

Report

General

Right

Left

Impression + Reco

Final Report

← Back

Patient Name: sofia

Patient ID: P100016

Gender: Female

DOB: August 4, 1997 (28 Y)

Patient Form: S Form

Date: 18 Sept 2025

Scan Center: MSC

Other reports

S.No	Report Title	File
No Reports Found		

Other Old reports

S.No	Report Title	File
1	Thermogram	View

Time stamp

Created by	Start	End	Status
sofia	18 Sept 2025	18 Sept 2025	Patient Intake For...

Final Report Preview

Ease QT 10.10 Auto Report ☒

Normal

B

I

U

NAME	sofia	DOB	August 4, 1997 (28 y)
GENDER	F	SCAN CENTER	MSC, Bangalore
USERID	P100016	DATE OF VISIT	Sep 18, 2025

QT ULTRASOUND BREAST IMAGING

Indication: Screening - Routine Annual Checkup.

History:

Reason for having this QT scan (Patient Stated): sss.

Patient is asymptomatic.

Imaging:

- sss.

TECHNIQUE: Transmission and reflection multiplanar 3-dimensional ultrasound imaging of both breasts was performed using the QT Ultrasound Series 2000 Model-A scanner.

Copyright © Wellthreen

Back

1/1

Next

z0SuYYZcbHUgC4gse+1qqr7Y0J7gUAEhCYITBAJM...

1 / 15

100%


1

S.No	Title	PAGE No
1	Introduction	2

Copyright © Wellthreen

The users can see uploaded Patient’s old reports and their current report.

ANALYTICS

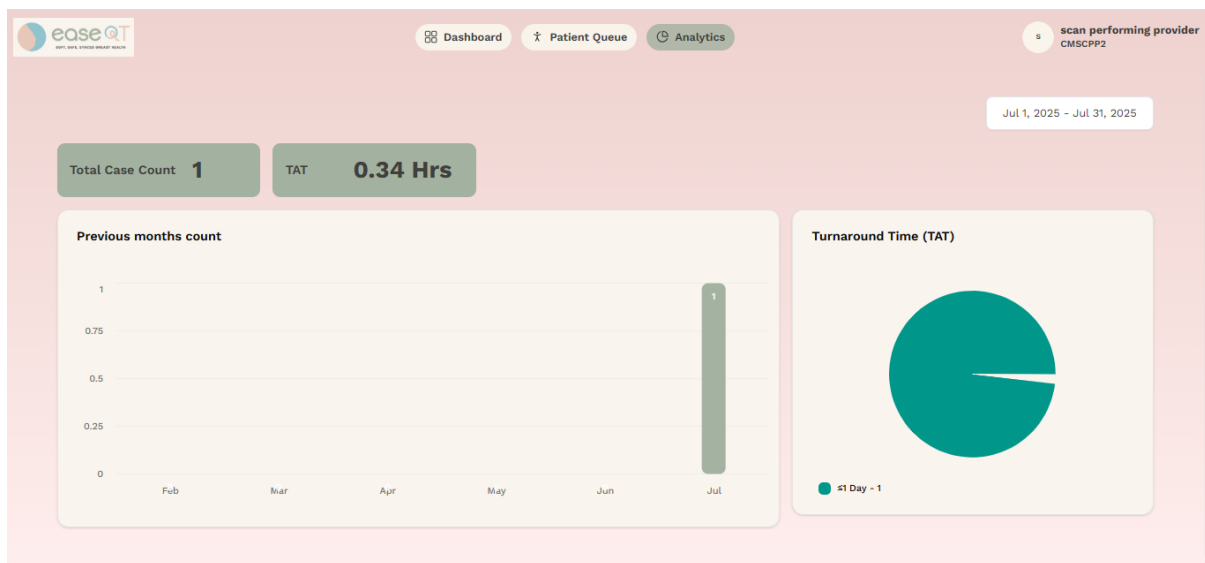

In analytics, the scan centre performing provider can view total case count, TAT, previous 6 months count.

GENERAL TIPS

1. If you open something in edit mode, save it before exiting-otherwise, others won't be able to use it.
2. After addressing pending remarks, leave another remark saying "Done" and set Assign to as None.
3. The Filter option is available for sorting.
4. For bulk downloads, use the checkboxes to select patients and then download in bulk.
5. To add a pending remark, click the + icon.
6. When viewing a report, always click Back to exit-otherwise, the case you opened may freeze.
7. Use Image Preview to add images to a report.
8. To refresh the page, press Ctrl + Shift + R.
9. For DICOM downloads and voice dictation, use Chrome or Edge. (For Firefox, click Save. Do not use Internet Explorer.)
10. Always log out properly-do not just close the browser or tab.