


ease @T


SOFT, SAFE, SYNCED BREAST HEALTH

RADIOLOGY USER MANUAL

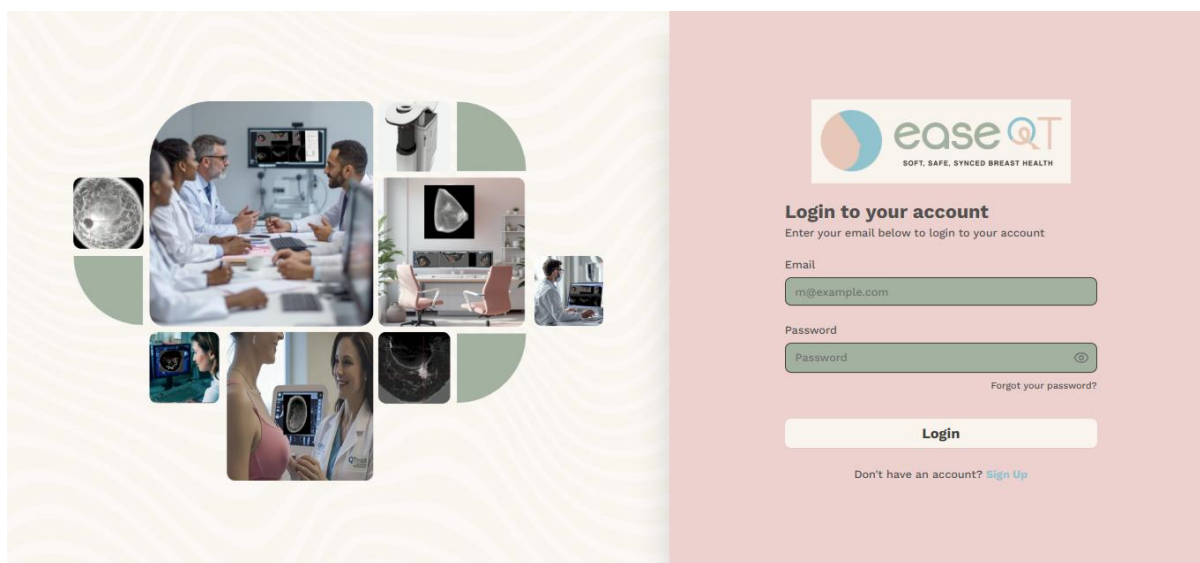

Once the scribe completes the report, the radiologist will begin their review, verify the scribe's entries, and proceed to finalize the report.

The radiology can sign in using the email and password provided by WellthGreen.

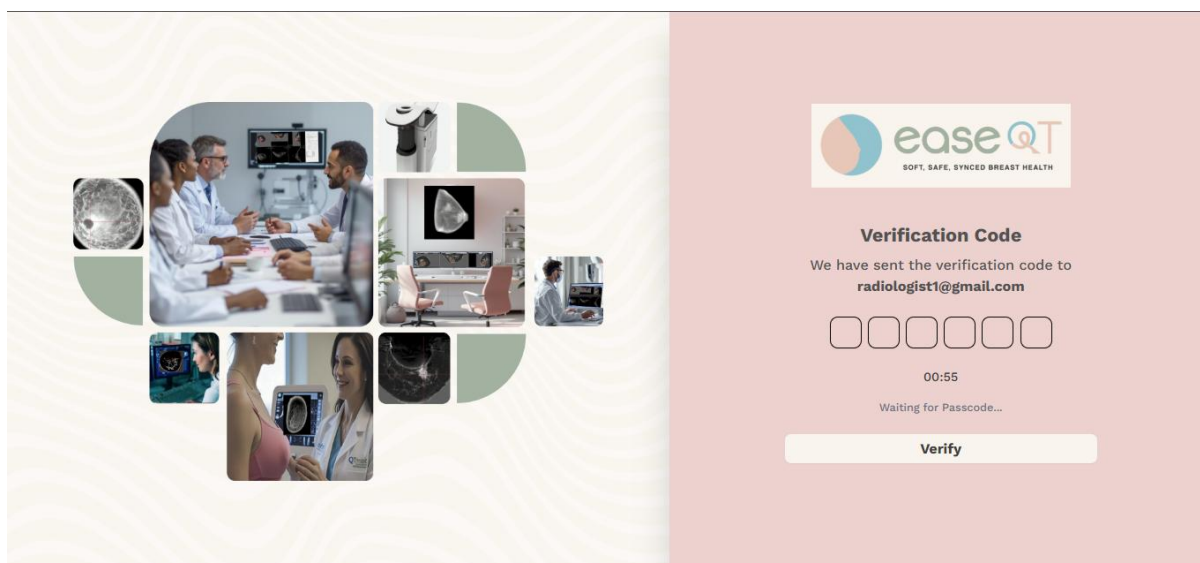

The Passcode will be generated for the verification.

Of Care. I understand that all such data will be handled securely and only by authorized personnel.

5. Credential Use and Responsibility

I confirm that all credentials, licenses, or certifications I provide during registration are accurate and current. I agree to keep my login credentials confidential and to notify the platform administrator of any suspected breach or unauthorized access.

6. Communication & Updates

I consent to receiving essential updates, notifications, and communications related to the EaseQT platform via email or in-app alerts.

☒ I agree to the above universal terms and policies.


7. Radiologist Consent


- I will review and input findings based on available imaging and history.
- I acknowledge that my reports may be reviewed or finalized by a signatory doctor.
- I will adhere to the platform's turnaround time (TAT) policy.
- I agree to collaborate respectfully with peers and escalate unclear cases as needed.

☒ I agree to the Co-Reporting Doctor role-specific terms.

[Submit](#)

The radiologist needs to sign the form to proceed.





Reset Your Password
E-Mail: radiologist1@gmail.com


New Password

- ☒ At least 8 characters
- ☒ Contains a number
- ☒ Special character (!@#%&^&*)


Confirm Password

[Reset Password](#)


The radiologist need to reset password for their security purpose.

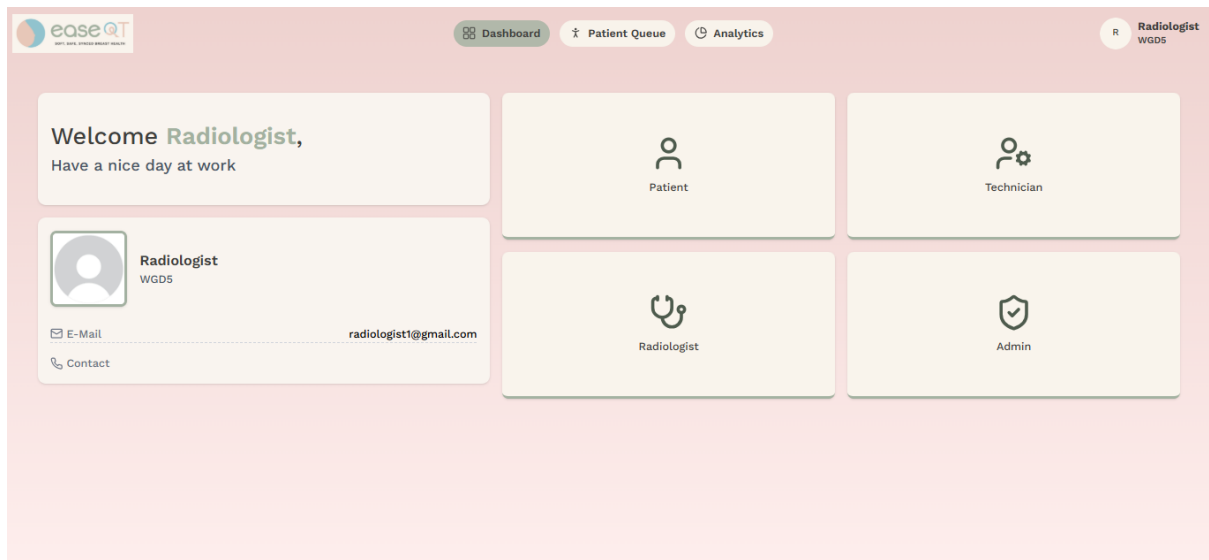


In the dashboard, the radiologist can view the patient brochure, technician guidelines, radiologist training material, and invoices.

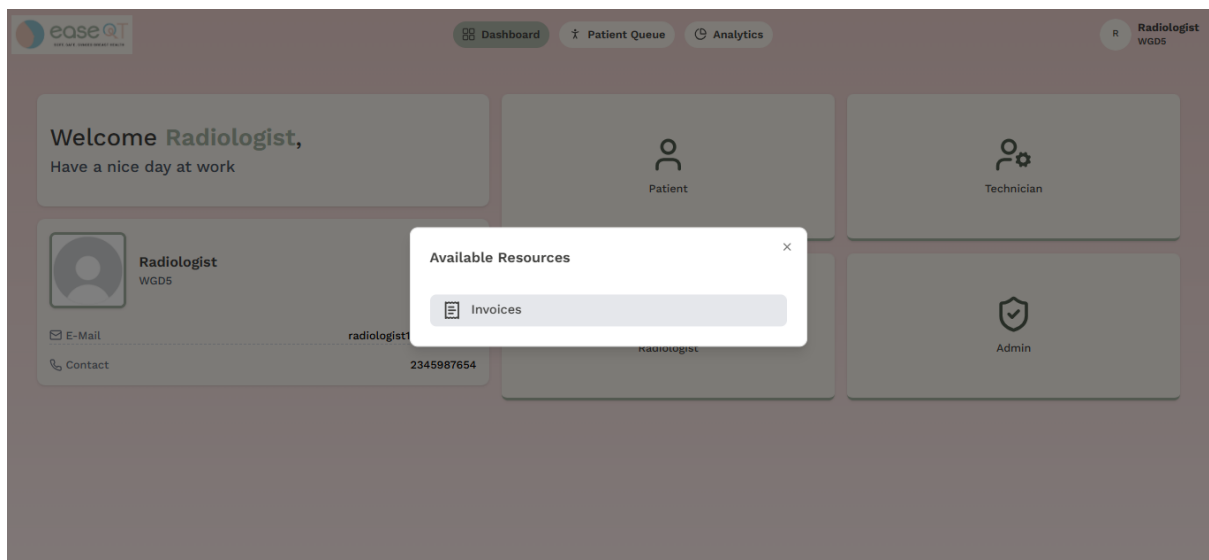


1. Invoices

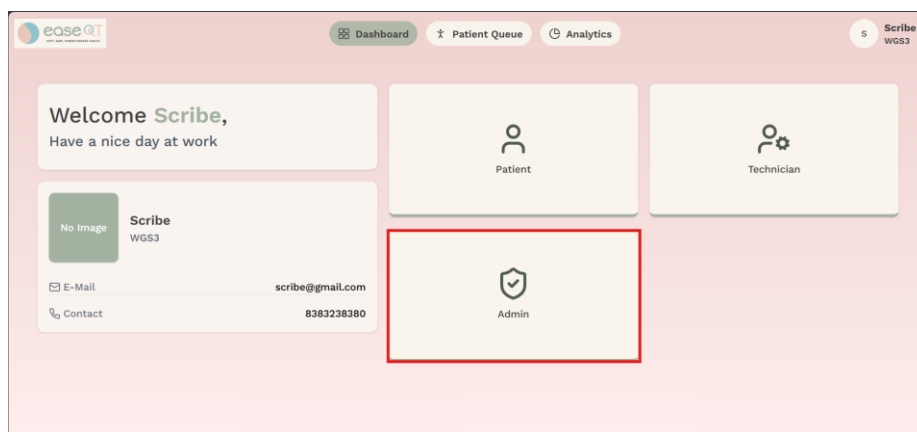


Clicking on the 'Admin' button triggers a popup window.

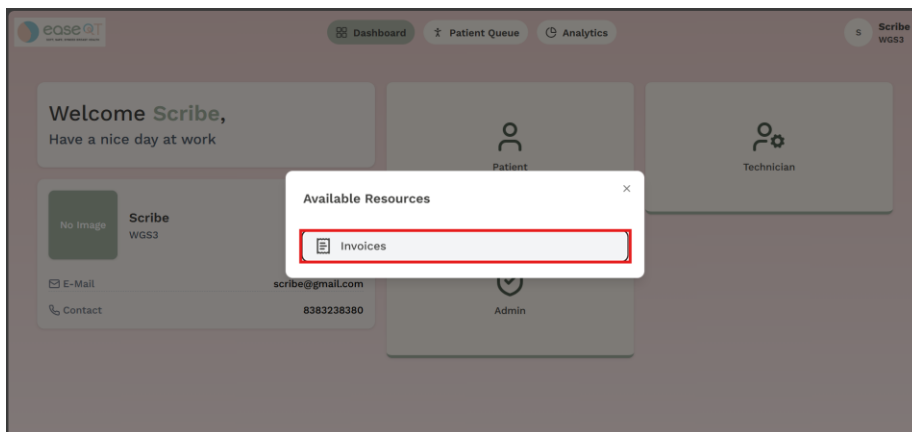

Clicking the 'Invoices' button open to the Invoice page.

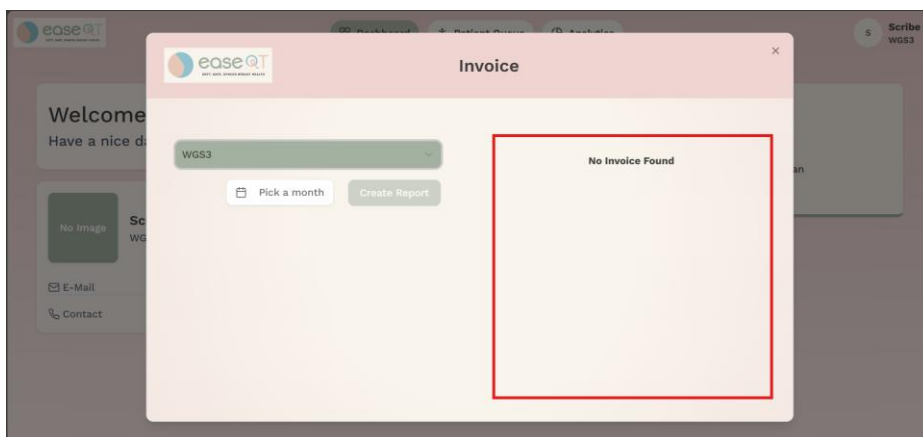

This section displays a list of invoices generated by you or the admin.

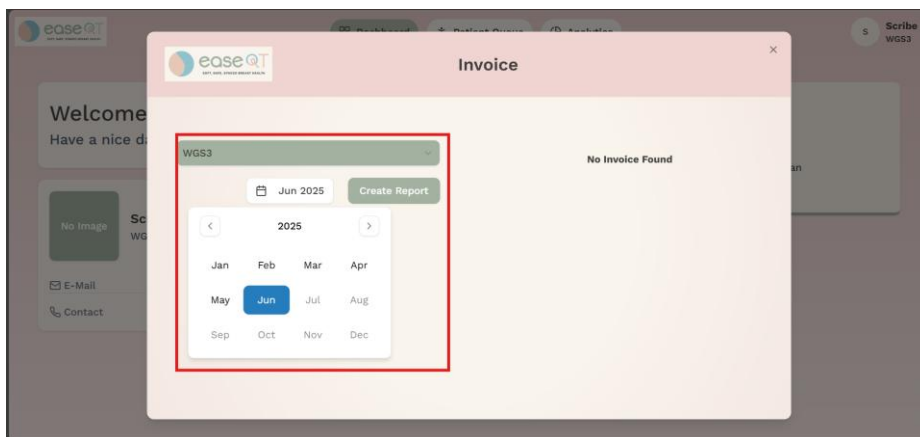

If you need an invoice, select the desired month from the enabled options and click the 'Create Report' button to proceed to the invoice creation page.

Invoice

← Back

From

Scribe
8382238380
scribe@gmail.com
PAN
GST
Address

Billing Period

June 1st, 2025
To
June 30th, 2025

Invoice Details

Quantity: 0
Amount (INR): 3000
Total: 0/-

To

Wellth Green Health Care Pvt. Ltd.
Address

UPI
UPI ID

Generate Report

The marked fields are auto-populated with data from the system and cannot be edited. Unmarked fields are editable. Once all the required details are filled in, click the 'Generate Report' button to create the invoice.

Invoice

WGS3
2025-06
Download

Pick a month
Create Report

Download

Generated invoices are listed here. If you need a copy, click the 'Download' button to download your invoice.

Invoice

No Invoice: 017013
Date: 2025-07-24
Billing Period: 2025-06-01 to 2025-06-30
Mode of Payment: BANK TRANSFER
Holder Name: Scribe
Number: 1234567890
Bank: Indian Bank
Branch: Salem
IFSC: #AIZ3456

From:
Name: Scribe
Phone: 8382238380
Email: scribe@gmail.com
PAN: 0151933304
GSTIN: 1234567890
Address: Salem

To:
Name: Wellth Green Health Care Pvt. Ltd.
Address: Chennai

S.No.	Description	Quantity	Amount (INR)
1.	Number of Scans Reported	0	3000/-

Total: 0/-

Declaration: I hereby declare that the information provided above is accurate, and the services listed have been duly completed during the billing period.

Signature:
Scribe
Date: 2025-07-24

Invoice Preview

PATIENT QUEUE

easeQT

HELP DATA PATIENTS HISTORY REPORTS

Dashboard

Patient Queue

Analytics

0

W

Wellthgreen Radiologist

WGBS

WellthGreen Radiologist

Download Dicom

Clear All Filters

Search all columns...

	Date	Scan Centre	Patient ID	Patient Name	Form	Tech Form	DICOM		DICOM	Old Report	QT Report	Report Status	Assign	Pending Remarks	Remarks
							L	R							
<input type="checkbox"/>	Sep 17, 2025	MSC	P100006	Testing12	Dc - View	View	(1)	-	View DICOM		View	SIGNED OFF (A)	<div><div></div></div>	No remarks...	<div><div>+</div><div></div></div>
<input type="checkbox"/>	Sep 18, 2025	MSC	P100006	Testing12	S - View	View	(1)	-	View DICOM		-	REVIEWED 2	<div>Assign <div></div></div>	R1 R2 stuck	<div><div>+</div><div></div></div>
<input type="checkbox"/>	Sep 22, 2025	MSC	P100006	Testing12	S - View	View	(1)	-	View DICOM		View	SIGNED OFF (A)	<div>Assign <div></div></div>	No remarks...	<div><div>+</div><div></div></div>

Rows per page

50

<<

< Previous

1

Next >

>>

Page 1 of 1

Copyright © Wellthgreen

The radiologist can view the patient form, technician form, and the DICOM files uploaded by the technician. Using these, they can start creating the report. They can also assign the report to another user if needed and add pending remark and can view previous remark also. Click on **View** to proceed with filling out the report. And they can multi select and download DICOM. They can upload and download the Patient's old report.

← Back

Patient Name: sofia
Patient ID: P100016
Gender: Female
DOB: August 4, 1997 (28 Y)
Patient Form: S Form
Date: 18 Sept 2025
Scan Center: MSC

Other reports

S.No	Report Title	File
No Reports Found		

Other Old reports

S.No	Report Title	File
No Reports Found		

Time stamp

Created by	Start	End	Status
sofia	18 Sept 2025	18 Sept 2025	Patient Intake Form

☒ Right Breast

Impression

0	Incomplete/Indeterminate Study
1	Likely Benign Findings
1a	Likely Benign Findings + Architectural Distortion
2	Low-Risk Characteristics
2a	Low risk - too small
3	Indeterminate Findings: Non palpable mass
3a	Indeterminate: Not Definitively Abnormal
3b	Indeterminate: Calcifications
3c	Indeterminate: Cysts
3d	Indeterminate : Postop changes + fluid/ seroma
3e	Indeterminate: Post-surgical vs. Recurrence
3f	Indeterminate: Likely Fibrosis
3g	Indeterminate: Long Standing Findings
4	Pain But No LUMP + QT Finding
4a	Palp Lump - No QT Finding or Ab/Normal image findings
4b	Palpable Lump + QT Finding

Recommendation

0	Repeat QT Imaging
1	Routine Annual Screening
1a	Routine Annual Screening
2	Short-Term Follow-Up 6m
2a	Short-term Follow-up 6m
3	Growth Rate Assessment 3m
3a	Confirm Benign Nature with Targ USG
3b	Targ USG for Cystic Mass Evaluation
3c	Clinical consultation with a breast specialist or health care provider for Calcification Management
3d	Clinical Consultation
3e	Clinical consultation + Surveillance
3f	Clinical Consultation + Surveillance
3g	If persists - Targ Usg
4	Clinical Monitoring
4a	Breast specialist Evaluation, neg never preclude/Targ USG.

The radiology can give the impression and recommendation by clicking on it. They can add extra impression and recommendation.

Time stamp

Created by	Start	End	Status
sofia	18 Sept 2025 16:45	18 Sept 2025 04:16	Patient Intake Form Filled
Technician	18 Sept 2025 04:17	18 Sept 2025 04:18	Technology Form Filled
Wellthgreen	18 Sept 2025 04:18	18 Sept 2025 04:18	Saved

Status saver

Load Template Reset to Default

Predraft Draft

Reviewed 1 Correct Reviewed 1 Edit

Reviewed 2 Correct Reviewed 2 Edit

Insert Signature Sign Off

Final Report Preview

Ease QT 10.10 Auto Report ☒

Normal **B I U**

The rest of the breast tissue appears unremarkable.

RIGHT BREAST:
 Impression:
 Recommendation:


LEFT BREAST:
 Impression:
 Recommendation:

Patients are encouraged to continue routine annual breast cancer screening as appropriate for their age and risk profile. In addition to imaging, monthly breast self-examinations are recommended. Patients should monitor for any new lumps, focal thickening, changes in breast size or shape, skin dimpling, nipple inversion or discharge, or any other unusual changes. If any new symptoms or palpable abnormalities are identified between scheduled screenings, patients should promptly consult their healthcare provider for further evaluation.

Copyright © Wellthgreen

In the report section, the radiologist will complete the report, with some fields auto-populated from the technician form. If they want to make corrections to the final report, they can switch from “Unsync (easeqt 10.10 option)” to “Sync(easeqt 10.10 option)” to enable editing. And they can submit by giving draft.

ANALYTICS

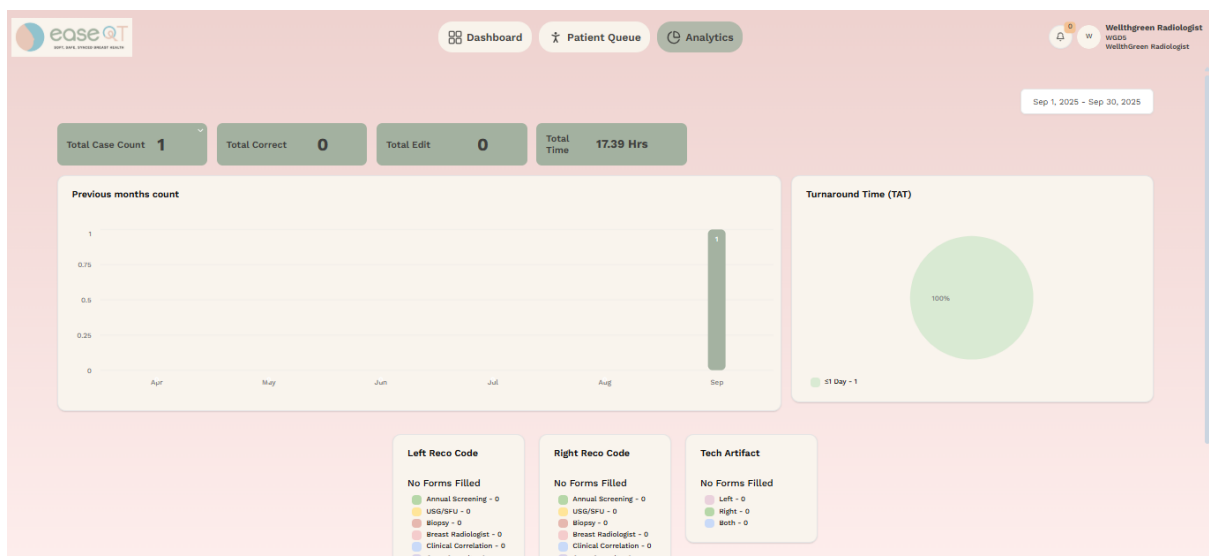

In Analytics, the radiologist can view their total case count, the number of correct and edited cases, TAT (Turnaround Time), data for the previous six months Artifacts and Reco code.