


ease @T


**SOFT, SAFE, SYNCED BREAST HEALTH**

## WELLGREEN PERFORMIN PROVIDER LOGIN USER MANUAL

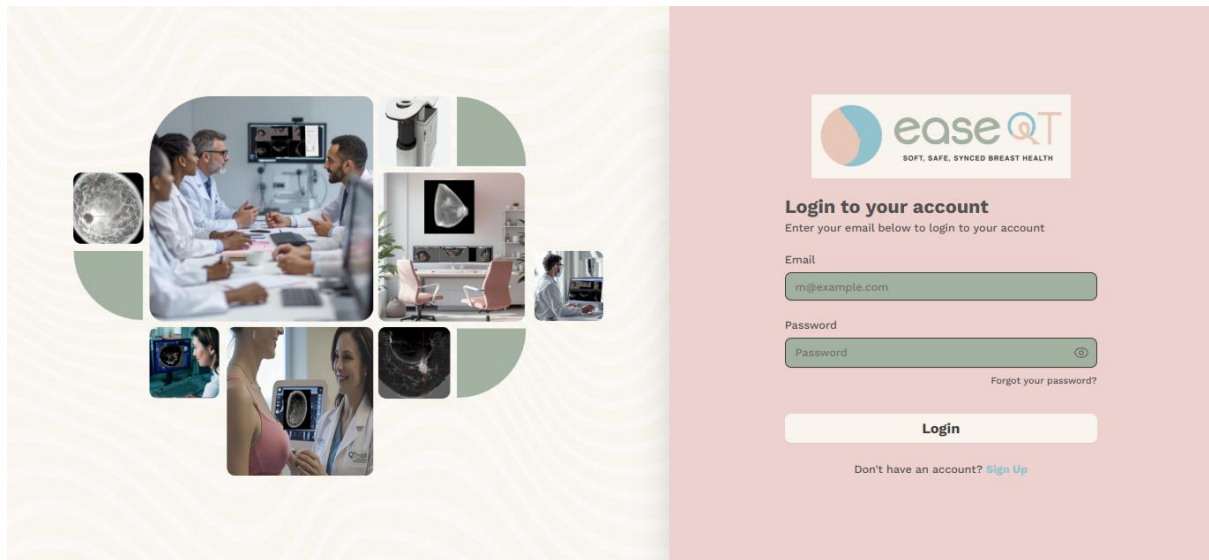

If master admin assigns to performing provider. They can start doing the report.

The performing provider can sign in using the email and password provided by WellthGreen.

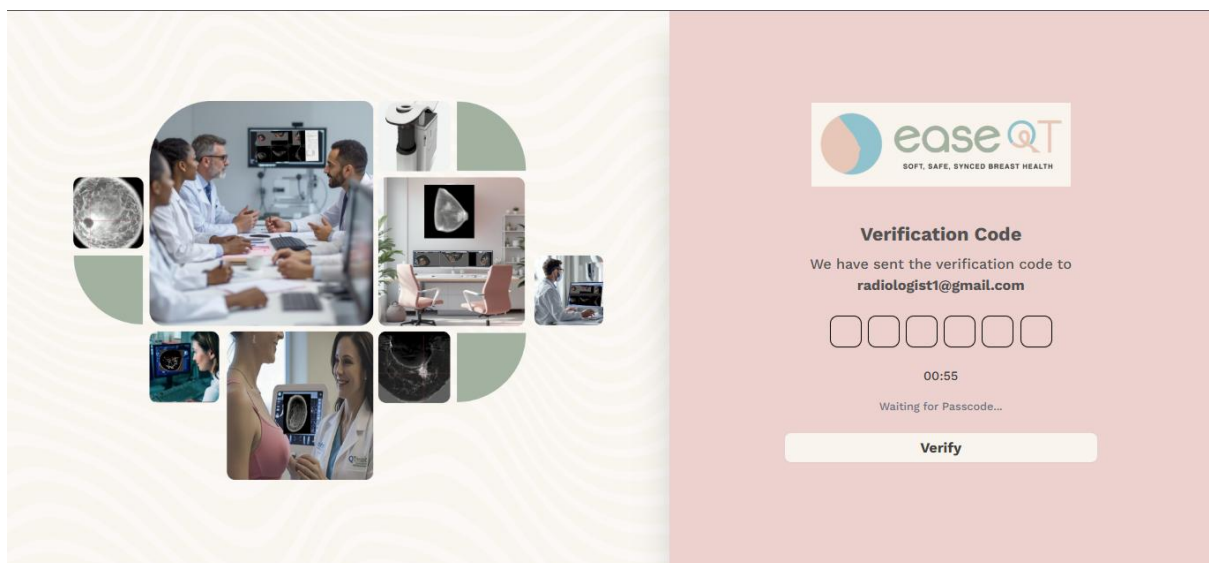


The verification code will be generated for the verification.



The image shows a 'User Consent Form' for the 'EaseQT Platform'. The form is presented in a window with a pink header. It contains four sections of consent text:

- 1. Consent to Platform Use**  
I hereby acknowledge that I am voluntarily registering to use the EaseQT platform for purposes directly related to clinical imaging, workflow management, or reporting.
- 2. Data Privacy & Security**  
I understand that my personal information and any patient-related data accessed through this platform may be securely stored and processed in compliance with global privacy laws (including HIPAA, DPA 2021, and GDPR where applicable). I agree not to share, download, or misuse any protected information obtained via the platform.
- 3. Confidentiality Agreement**  
I agree to maintain full confidentiality of any data I access, including but not limited to patient records, reports, clinical communications, and platform operations. Any breach of confidentiality may result in legal and professional consequences.
- 4. Data Consent**  
I consent to the transfer, storage, and processing of data as necessary for platform function and continuity of care. I understand that all such data will be handled securely and only by authorized personnel.


The performing provider needs to sign the form to proceed.

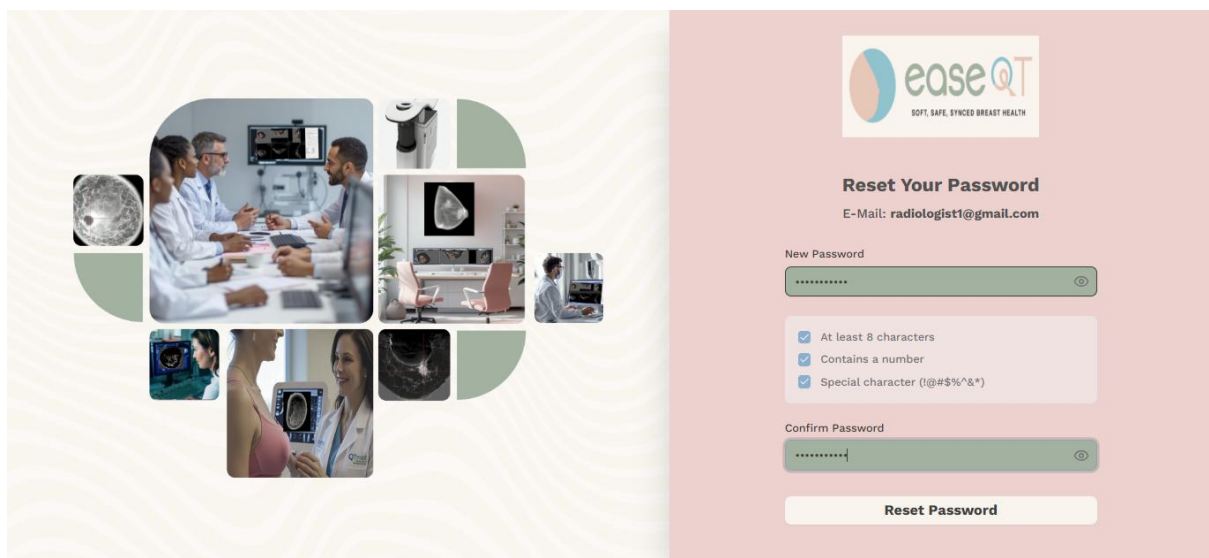


The image shows a 'Reset Your Password' form for the 'EaseQT' platform. The form is presented in a window with a pink header. It includes the EaseQT logo and the tagline 'SOFT, SAFE, SYNCED BREAST HEALTH'. The form contains the following fields and elements:


- Reset Your Password**
- E-Mail: radiologist1@gmail.com
- New Password field (password type)
- Confirm Password field (password type)
- Reset Password button
- Validation requirements (all checked):
  - At least 8 characters
  - Contains a number
  - Special character (!@#%&^&\*)

The performing provider needs to reset password for their security purpose.

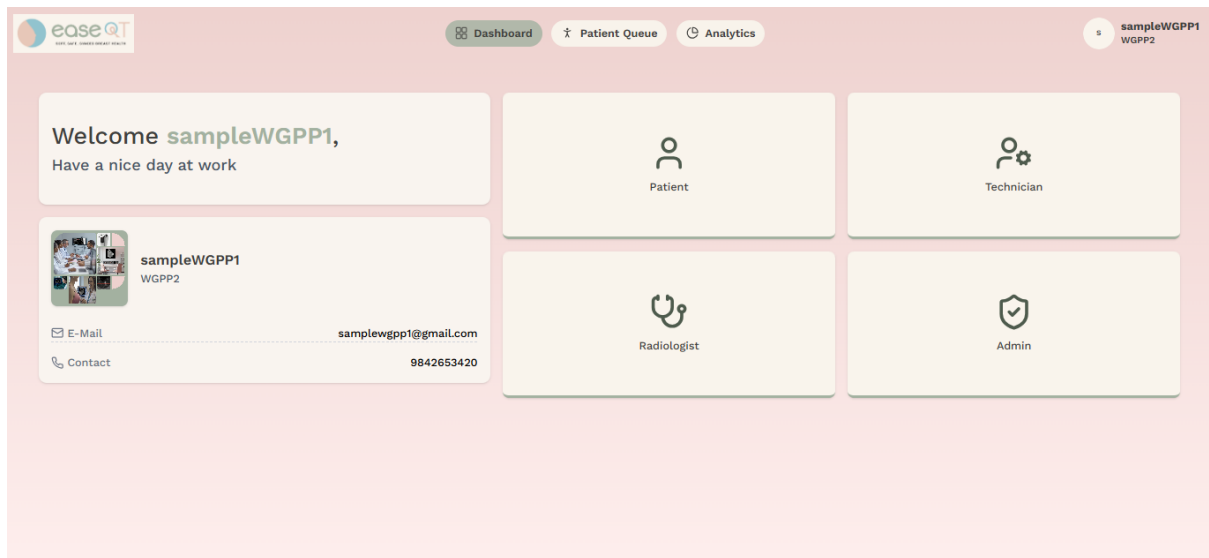


In the dashboard, the performing provider can view the patient brochure, general guidelines, disclaimer, patient consent form, technician guidelines and technician consent form. They can also view the radiologist training material, radiologist consent form and invoices.

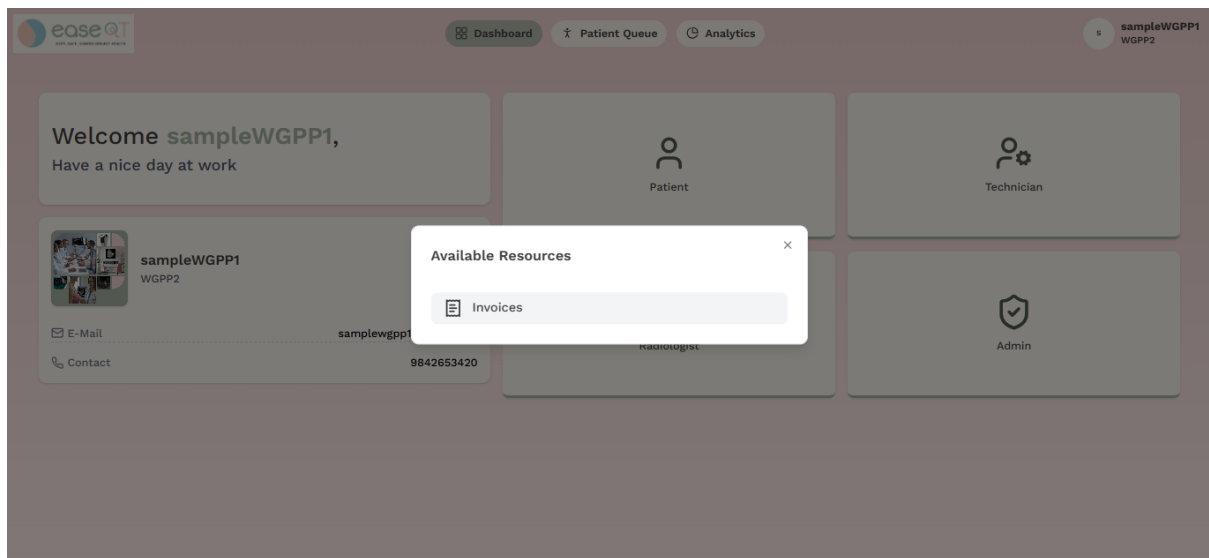

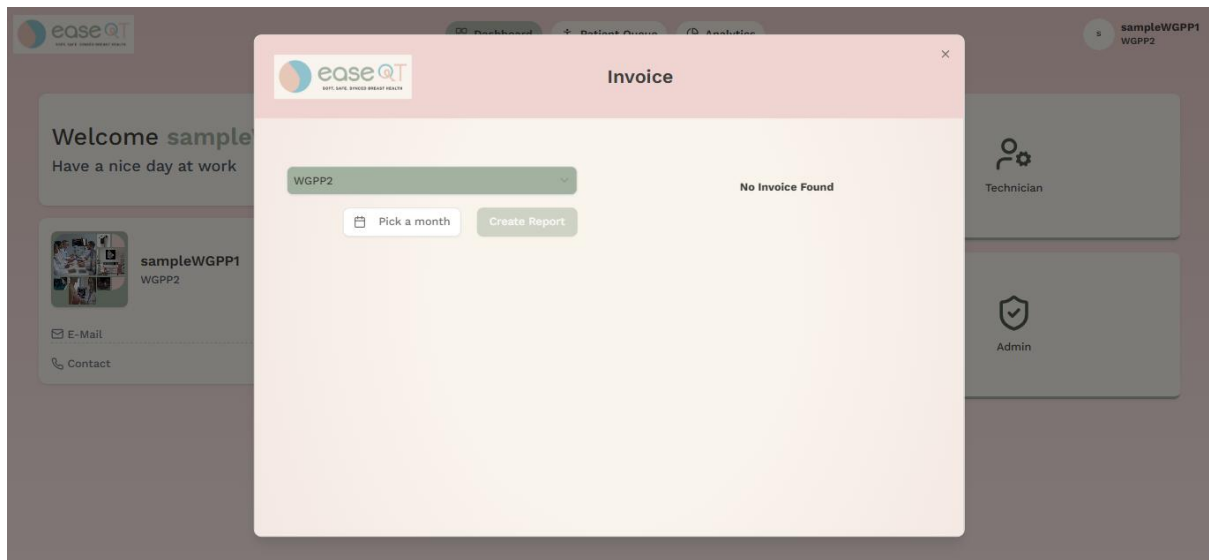


The performing provider will create their own invoice by clicking on admin. they need click on which month and create report and give necessary detail and generate report


## 1. Invoices

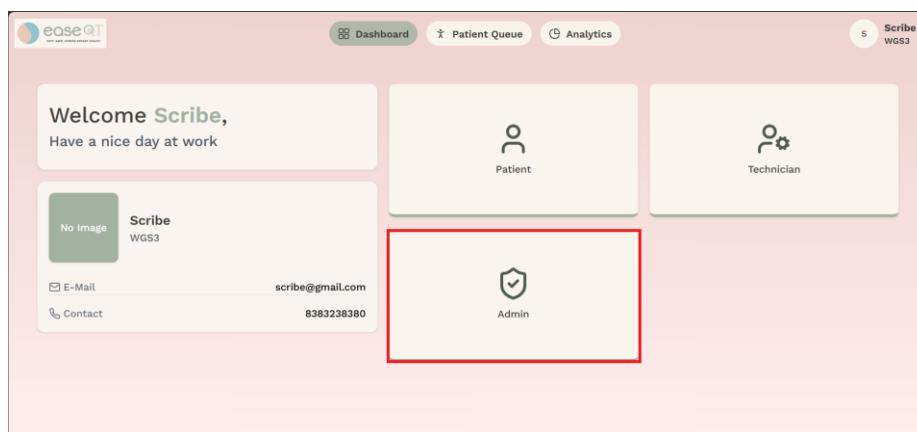


Clicking on the 'Admin' button triggers a popup window.

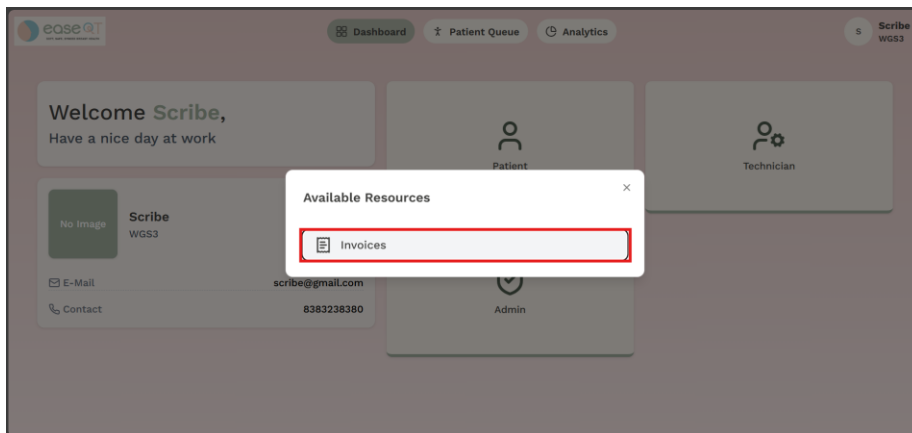

Clicking the 'Invoices' button open to the Invoice page.

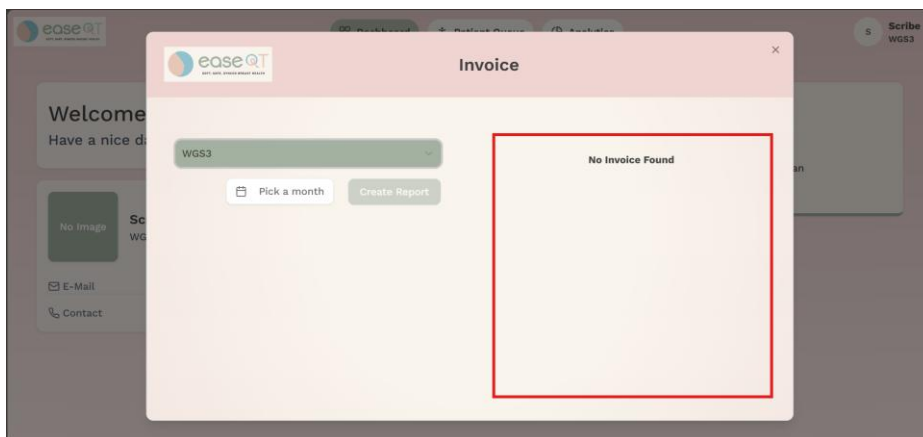

This section displays a list of invoices generated by you or the admin.

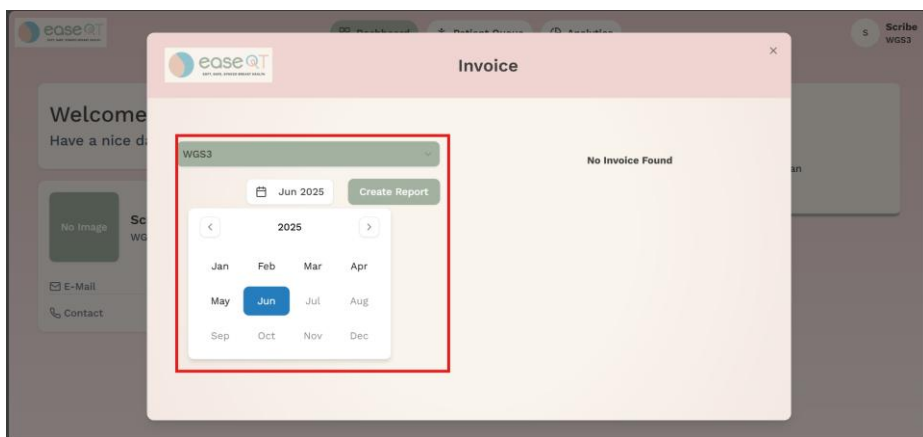

If you need an invoice, select the desired month from the enabled options and click the 'Create Report' button to proceed to the invoice creation page.

**Invoice**

← Back

**From**

Scribe  
8382238380  
scribe@gmail.com  
PAN  
GST  
Address

**Billing Period**

June 1st, 2025  
To  
June 30th, 2025

**Invoice Details**

Quantity: 0  
Amount (INR): 3000  
Total: 0/-

**To**

Wellth Green Health Care Pvt. Ltd.  
Address

UPI  
UPI ID

Generate Report

The marked fields are auto-populated with data from the system and cannot be edited. Unmarked fields are editable. Once all the required details are filled in, click the 'Generate Report' button to create the invoice.

**Invoice**

WGS3

2025-06

Pick a month

Create Report

Generated invoices are listed here. If you need a copy, click the 'Download' button to download your invoice.

**Invoice**

No Invoice: 07003  
Date: 2025-07-24  
Billing Period: 2025-06-01 to 2025-06-30  
Mode of Payment: BANK TRANSFER  
Holder Name: Scribe  
Number: 1234567890  
Bank: Indian Bank  
Branch: Salem  
IFSC: IAI23456

**From:**  
Name: Scribe  
Phone: 8382238380  
Email: scribe@gmail.com  
PAN: DT51933304  
GSTIN: 1234567890  
Address: Salem

**To:**  
Name: Wellth Green Health Care Pvt. Ltd.  
Address: Chennai

S.No.	Description	Quantity	Amount (INR)
1.	Number of Scans Reported	0	3000/-


**Total: 0/-**

Declaration: I hereby declare that the information provided above is accurate, and the services listed have been duly completed during the billing period.

Signature  
Scribe  
Date: 2025-07-24

Invoice Preview

## PATIENT QUEUE



Dashboard

Patient Queue

Analytics

0

W

Wellthgreen Performing Provider  
WGPP2  
WellthGreen Performing Provider

Download Dicom

Download Final Reports

Download Patient Consent

Clear All Filters

Search all columns...

<input type="checkbox"/>	Date	Scan Centre	Patient ID	Patient Name	Consent	Form	Tech Form	DICOM		Old Report	QT Report	Report Status	Assign	Pending Remarks	Remarks
								L	R						
<input type="checkbox"/>	Sep 17, 2025	MSC	P100006	Testing12	View	Dc - View	View	(1)	-		View	SIGNED OFF (A)	AD10000	No remarks...	<div>+ </div>
<input type="checkbox"/>	Sep 18, 2025	MSC	P100006	Testing12	View	S - View	View	(1)	-		View	REVIEWED 2	Assign	R1 R2 stuck	<div>+ </div>
<input type="checkbox"/>	Sep 22, 2025	MSC	P100006	Testing12	View	S - View	View	(1)	-		View	SIGNED OFF (A)	Assign	No remarks...	<div>+ </div>
<input type="checkbox"/>	Sep 18, 2025	MSC	P100016	sofia	View	S - View	View	(1)	-		View	DRAFT	Assign	No remarks...	<div>+ </div>

Rows per page50

<<

< Previous

1

Next >

>>

Page 1 of 1

In the patient queue, the performing provider can view the patient list, patient forms, technician view and DICOM files. They can fill out reports, check the report status, assign the case to another user and add pending remarks and can view previous remark also. Click on **View** to proceed with filling out the report. And they can multi select and download DICOM and final reports. They can upload and download Patient's old report.

Created by	Start	End	Status
sofia	18 Sept 2025 16:45	18 Sept 2025 04:16	Patient Intake Form Filled
Technician	18 Sept 2025 04:17	18 Sept 2025 04:18	Technician Form Filled
Wellthgreen	18 Sept 2025 04:18	18 Sept 2025 04:18	Saved

### Final Report Preview

Ease QT 10.10 Auto Report ☒

NAME	sofia	DOB	August 4, 1997 (28 y)
GENDER	F	SCAN CENTER	MSC, Bangalore
USERID	P100016	DATE OF VISIT	Sep 18, 2025

**QT ULTRASOUND BREAST IMAGING**  
 Indication: Screening : Routine Annual Checkup.  
 History:  
 Reason for having this QT scan (Patient Stated): sss.  
 Patient is asymptomatic.  
 Imaging:  
 • sss.

**TECHNIQUE:** Transmission and reflection multiplanar 3-dimensional ultrasound imaging of both breasts was performed using the QT Ultrasound Series 2000 Model-A scanner.

**Status saver**

Load Template Reset to Default

Predraft Draft


Reviewed 1 Correct Reviewed 1 Edit

Reviewed 2 Correct Reviewed 2 Edit

Insert Signature Sign Off

In the report section, the performing provider can review and make changes to reports completed by the radiologist. There are “Reviewed 1 Correct” and “Reviewed 2 Edit” buttons available for the performing provider. If the radiologist has completed the report correctly, the performing provider will select “Reviewed 1 Correct.” If it is incorrect, the performing provider will select “Reviewed 2 Edit.” After verifying the report, the performing provider





EASE QT

SAFE, SIMPLE, SENSITIVE BREAST HEALTH

Patient Form

Technician Form

Dicom

Report

General

Right

Left

Impression + Reco

Final Report

← Back

Patient Name: sofia

Patient ID: P100016

Gender: Female

DOB: August 4, 1997 (28 Y)

Patient Form: S Form

Date: 18 Sept 2025

Scan Center: MSC

Other reports

S.No	Report Title	File
No Reports Found		

Other Old reports

S.No	Report Title	File
1	Thermogram	<a href="#">View</a>

Time stamp

Created by	Start	End	Status
sofia	18 Sept 2025	18 Sept 2025	Patient Intake For

Final Report Preview

Ease QT 10.10 Auto Report ☒

Normal

↺

↻

B

I

U

≡

≡


≡

≡

A

A

📷



NAME	sofia	DOB	August 4, 1997 (28 y)
GENDER	F	SCAN CENTER	MSC, Bangalore
USERID	P100016	DATE OF VISIT	Sep 18, 2025

QT ULTRASOUND BREAST IMAGING

Indication: Screening : Routine Annual Checkup.

History:


Reason for having this QT scan (Patient Stated): sss.

Patient is asymptomatic.

Imaging:

\* sss.

TECHNIQUE: Transmission and reflection multiplanar 3-dimensional ultrasound imaging of both breasts was performed using the QT Ultrasound Series 2000 Model-A scanner.

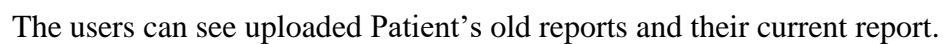

Patient Form

Technician Form

Dicom

Report

General

Right

Left

Impression + Reco

Final Report

← Back

Patient Name: sofia

Patient ID: P100016

Gender: Female

DOB: August 4, 1997 (28 Y)

Patient Form: S Form

Date: 18 Sept 2025

Scan Center: MSC

Other reports

S.No	Report Title	File
No Reports Found		

Other Old reports

S.No	Report Title	File
1	Thermogram	<a href="#">View</a>
2	Mammogram	<a href="#">View</a>
3	Breast MRI	<a href="#">View</a>

Time stamp

Created by	Start	End	Status
------------	-------	-----	--------

Right Breast

Impression

0	Incomplete/Indeterminate Study
1	Likely Benign Findings
1a	Likely Benign Findings + Architectural Distortion
2	Low-Risk Characteristics
2a	Low risk - too small
3	Indeterminate Findings: Non palpable mass
3a	Indeterminate: Not Definitively Abnormal
3b	Indeterminate: Calcifications
3c	Indeterminate: Cysts
3d	Indeterminate : Postop changes + fluid/ seroma
3e	Indeterminate: Post-surgical vs. Recurrence
3f	Indeterminate: Likely Fibrosis
3g	Indeterminate: Long Standing Findings
4	Pain But No LUMP + QT Finding
4a	Palp Lump - No QT Finding or Ab/Normal image findings
4b	Palpable Lump + QT Finding

Recommendation

0	Repeat QT Imaging
1	Routine Annual Screening
1a	Routine Annual Screening
2	Short-Term Follow-Up 6m
2a	Short-term Follow-up 6m
3	Growth Rate Assessment 3m
3a	Confirm Benign Nature with Targ USG
3b	Targ USG for Cystic Mass Evaluation
3c	Clinical consultation with a breast specialist or health care provider for Calcification Management
3d	Clinical Consultation
3e	Clinical consultation + Surveillance
3f	Clinical Consultation + Surveillance
3g	If persists - Targ Usg
4	Clinical Monitoring
4a	Breat specialist Evaluation, neg never preclude/Targ USG.

Copyright © Wellthgreen

The performing provider can give the impression and recommendation by clicking on it. They can add extra impression and recommendation.

Patient Form

Technician Form

Dicom

Report

General

Right

Left

Impression

Final Report

Created by	Start	End	Status
Admin	Aug 2025 03:14	Aug 2025 02:35	Changes Completed
Admin	11 Aug 2025 02:35	11 Aug 2025 03:13	Changes Completed
testing performing provider	11 Aug 2025 03:13	-	-

Load Template

Reset to Default

Predraft

Draft

Reviewed 1 Correct

Reviewed 1 Edit

Reviewed 2 Correct

Reviewed 2 Edit

Insert Signature

Sign Off

Final Report Preview

Unsync

Normal

B I U

NAME

Alpha Testing Patient

AGE/GENDER

30 / F

USERID

P100015

DOB

2001-01-02

SCAN CENTER

MSC, nrxjkdvn

DATE OF VISIT

2026-01-03

QT ULTRASOUND BREAST IMAGING

HISTORY :

A 30 year old woman with IBIS Tyrer-Cuzick score of % and AURIA breast cancer test with mastectomy, lumpectomy, breast reconstruction, augmentation surgery done.

Biopsy - Result: Normal

Symptoms: Lump or thickening in left at nipple, 7'o Clock and right at nipple, 9'o Clock, 12'o Clock: Pain in left at nipple, 5'o Clock, 7'o Clock and right at nipple, 7'o Clock: Nipple changes in right/inverted: Nipple discharge in left and right

Thermogram taken on 2025-01-08 shows abnormal results. The report is not available.

Mammogram taken on 2025-08-04 shows normal results. The report is available.

Breast ultrasound / HERscan taken on 2025-07-16 shows normal results. The report is not available.


Breast MRI shows unknown results. The report is not available.

PET/CT Scan shows abnormal results. The report is not available.

QT Imaging taken on 2025-08-04 shows abnormal results. The report is available.

Copyright © Wellthgreen

There is an option called "Reset to Default." If the report was done incorrectly, they can use this option to revert it to the default state.



Patient Form

Technician Form

Dicom

Report

General

Right

Left

Impression

Final Report

← Back

Name: Alpha Testing Patient  
Age & Gender: 30 / Female  
Date: 03 Jan 2026  
Scan Center: MSC

Created by	Start	End	Status
Admin	07 Aug 2025 03:14	11 Aug 2025 02:35	Saved Changes
Admin	11 Aug 2025 02:35	11 Aug 2025 03:13	Saved Changes
testing performing provider	11 Aug 2025 -	-	-

Load Template

Reset to Default

Predraft

Draft


Reviewed 1

Final Report Preview

Ease QT 10.10 Auto Report

☒

Normal





NAME	Alpha Testing Patient	DOB	2001-01-02
AGE/GENDER	30 / F	SCAN CENTER	MSC, nmxjkcdvn
USERID	P100015	DATE OF VISIT	2026-01-03

### QT ULTRASOUND BREAST IMAGING

**HISTORY :**  
A 30 year old woman who underwent mastectomy, lumpectomy, breast reconstruction, augmentation surgery.  
Biopsy - Result: Normal  
Symptoms: Lump or thickening in left at nipple, 7'o Clock and right at nipple, 9'o Clock, 12'o Clock: Pain in left at nipple, 5'o Clock, 7'o Clock and right at nipple, 7'o Clock: Nipple changes in right inverted: Nipple discharge in left and right  
Thermogram taken on 2025-01-08 shows abnormal results. The report is not available.  
Mammogram taken on 2025-08-04 shows normal results. The report is available.  
Breast ultrasound / HERscan taken on 2025-07-16 shows normal results. The report is not available.  
Breast MRI shows unknown results. The report is not available.  
PET/CT Scan shows abnormal results. The report is not available.  
QT Imaging taken on 2025-08-04 shows abnormal results. The report is available.  
Other imaging of scans (like bone scans, Coonnect, Wellthreen) shows abnormal results. The report is not available.

The Wellthgreen Performing provider can view the full report here. If they want to make changes or load a template, they want sink that need to be on. If they want to access the report, they need to switch on the 10.10 auto-report feature.

## ANALYTICS

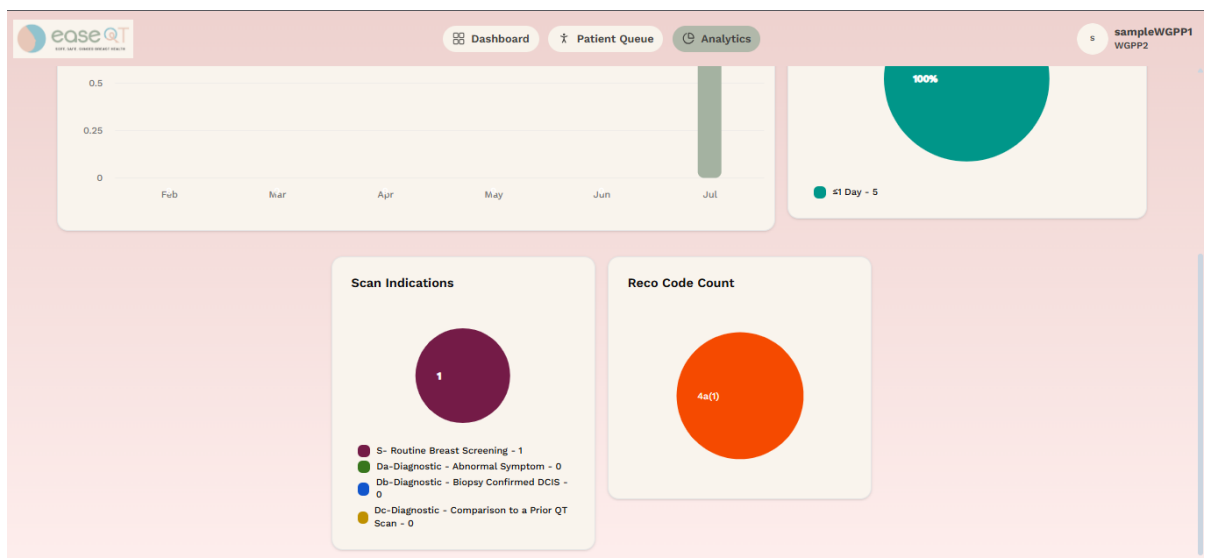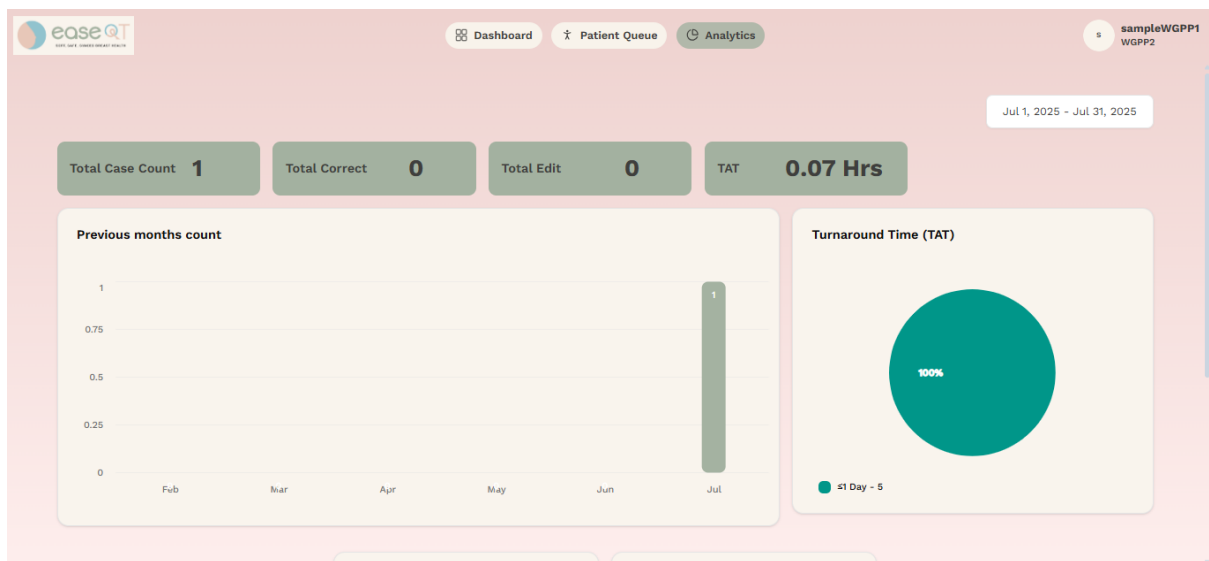

The admin can track analytics for all scan centres, users and their own activities. For their personal analytics, they can view the total count, number of edits, number of correct reports, counts for the previous six months, TAT (Turnaround Time), scan indications and recommendation codes.

## **GENERAL TIPS**

1. If you open something in edit mode, save it before exiting-otherwise, others won't be able to use it.
2. After addressing pending remarks, leave another remark saying "Done" and set Assign to as None.
3. The Filter option is available for sorting.
4. For bulk downloads, use the checkboxes to select patients and then download in bulk.
5. To add a pending remark, click the + icon.
6. When viewing a report, always click Back to exit-otherwise, the case you opened may freeze.
7. Use Image Preview to add images to a report.
8. To refresh the page, press Ctrl + Shift + R.
9. For DICOM downloads and voice dictation, use Chrome or Edge. (For Firefox, click Save. Do not use Internet Explorer.)
10. Always log out properly-do not just close the browser or tab.