





ease @T

SOFT, SAFE, SYNCED BREAST HEALTH


ADMIN USER MANUAL

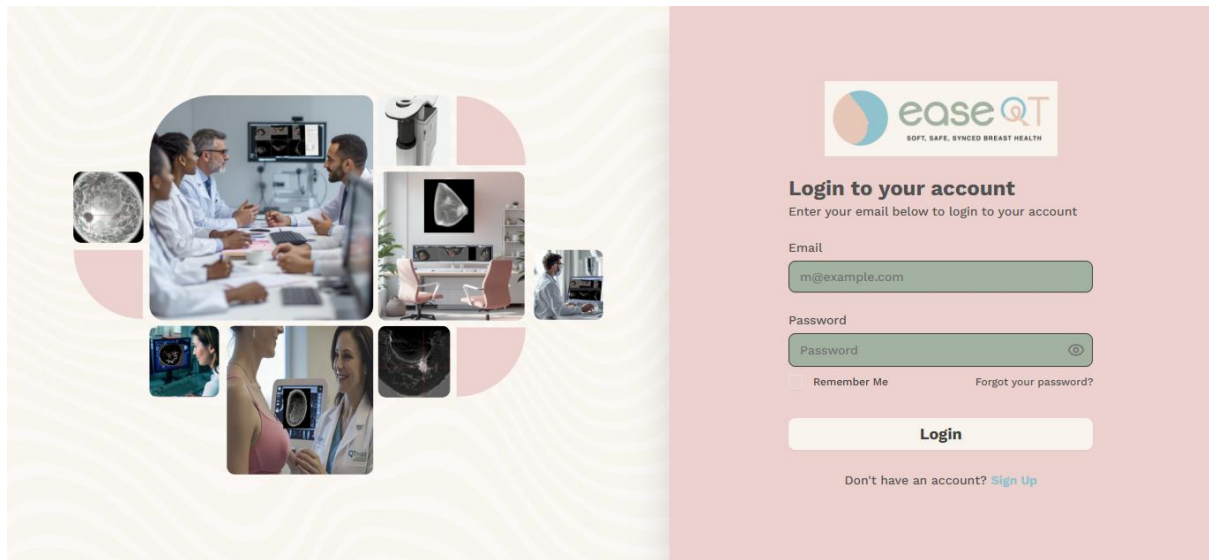

The admin can log in using the email and password provided by the ZAdroit team.

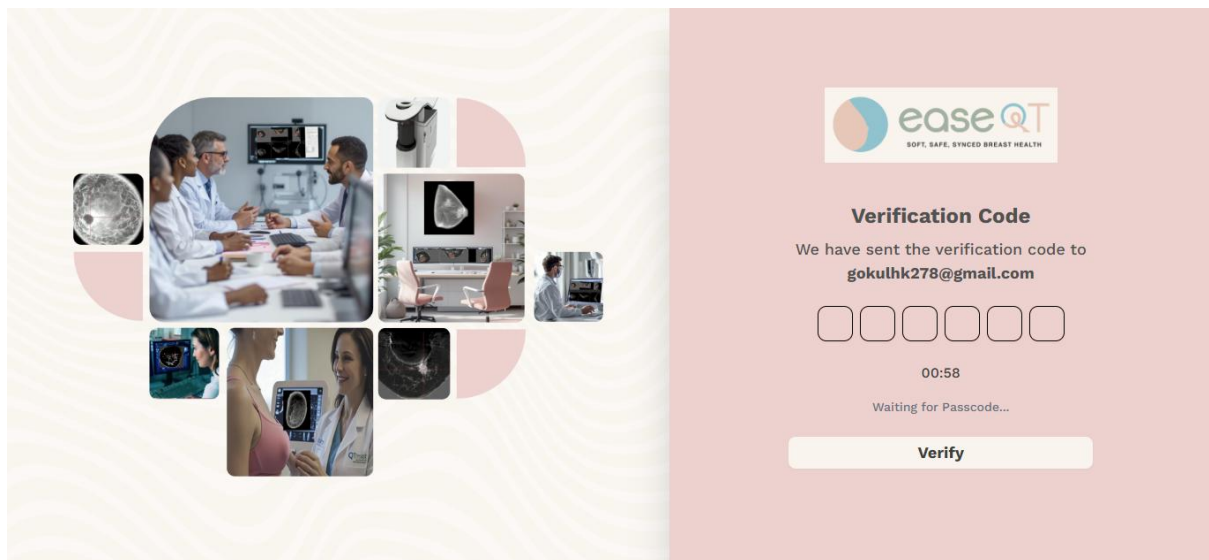


A verification code will be sent to the email for verification, and it needs to be entered in the application.

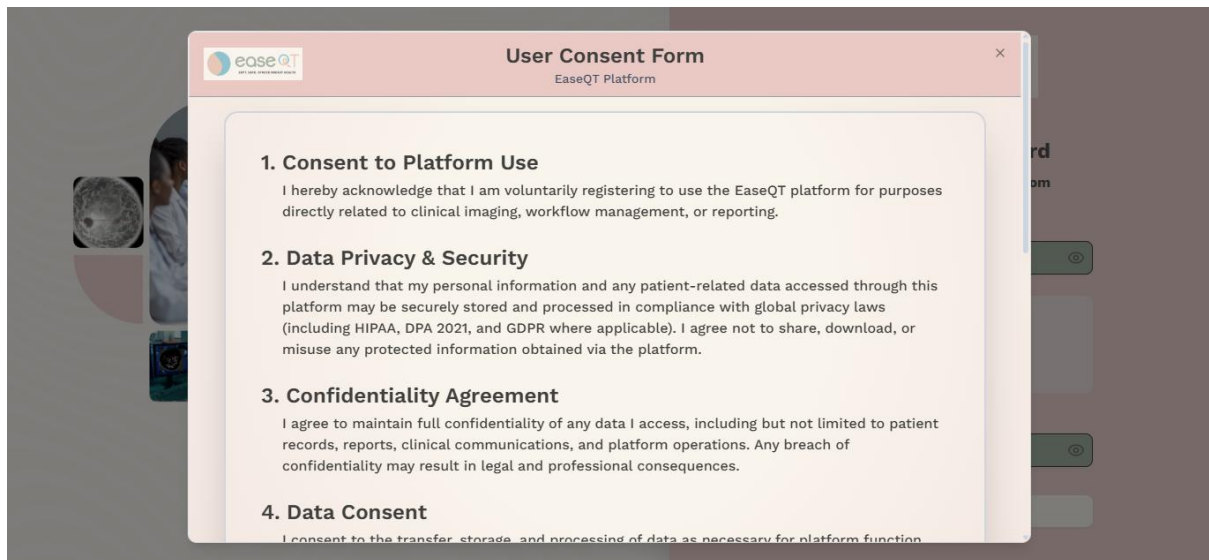

The image shows a 'User Consent Form' window for the 'EaseQT Platform'. The window has a pink header with the 'easeQT' logo and a close button. The form is divided into four sections:

- 1. Consent to Platform Use**
I hereby acknowledge that I am voluntarily registering to use the EaseQT platform for purposes directly related to clinical imaging, workflow management, or reporting.
- 2. Data Privacy & Security**
I understand that my personal information and any patient-related data accessed through this platform may be securely stored and processed in compliance with global privacy laws (including HIPAA, DPA 2021, and GDPR where applicable). I agree not to share, download, or misuse any protected information obtained via the platform.
- 3. Confidentiality Agreement**
I agree to maintain full confidentiality of any data I access, including but not limited to patient records, reports, clinical communications, and platform operations. Any breach of confidentiality may result in legal and professional consequences.
- 4. Data Consent**
I consent to the transfer, storage, and processing of data as necessary for platform function.

The admin needs to sign the form.

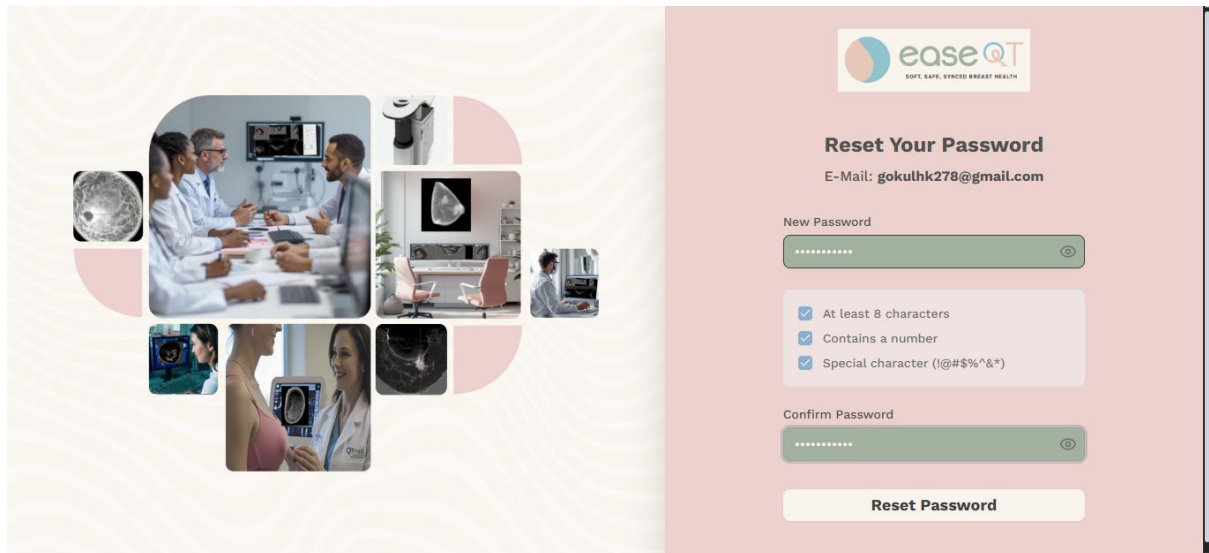


The image shows a 'Reset Your Password' form. On the left is a collage of medical-related images. The form itself is on a pink background and includes the 'easeQT' logo and the tagline 'SOFT. SAFE. SYNCHRONIZED BREAST HEALTH'. The title is 'Reset Your Password'. Below the title, the email 'E-Mail: gokulhk278@gmail.com' is displayed. The form has two password input fields: 'New Password' and 'Confirm Password', both with a strength indicator. Between the fields are three checkboxes, all of which are checked:

- ☒ At least 8 characters
- ☒ Contains a number
- ☒ Special character (!@#%&*)


At the bottom is a 'Reset Password' button.


The admin needs to reset their password.

DASHBOARD

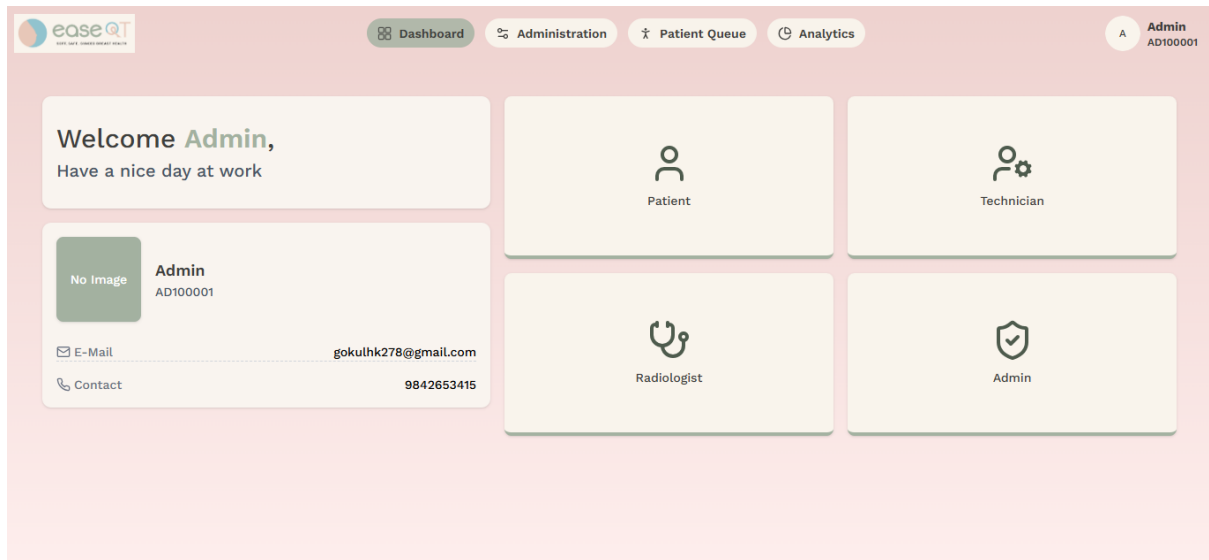


In the dashboard, the admin can view the patient brochure, general guidelines, disclaimer, patient consent form, technician guidelines and technician consent form. They can also view the radiologist training material, radiologist consent form and invoices.

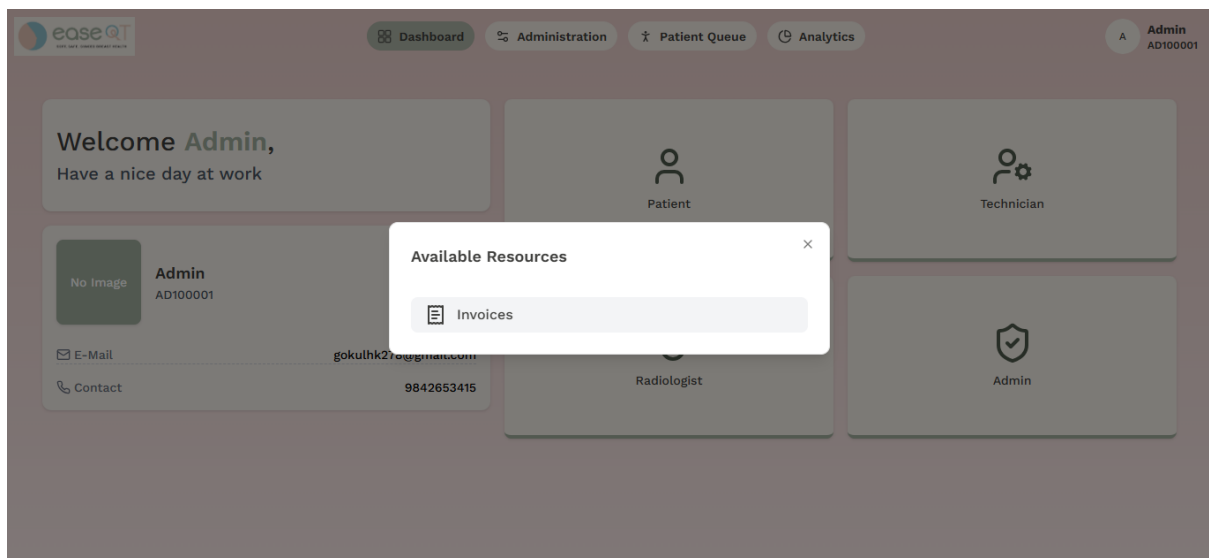

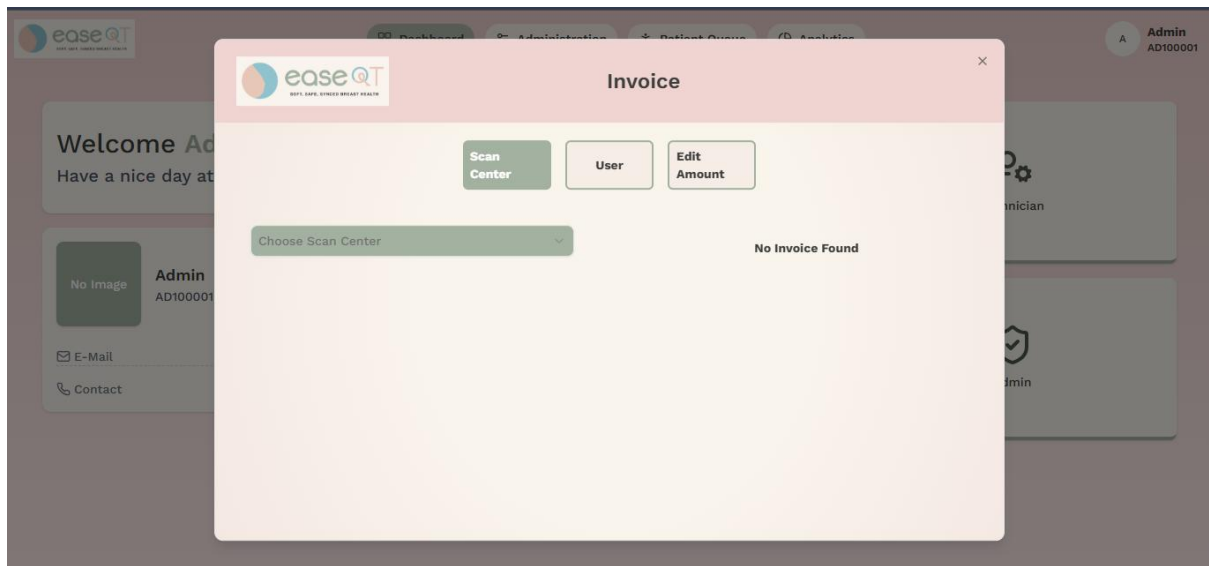


For invoice generation, the admin can select either a scan centre or a user, choose the desired month and click “Create Report.” A form will appear where the required details need to be entered, after which the report will be generated.

The marked fields are auto-populated with data from the system and cannot be edited. Unmarked fields are editable. Once all the required details are filled in, click the 'Generate Report' button to create the invoice.

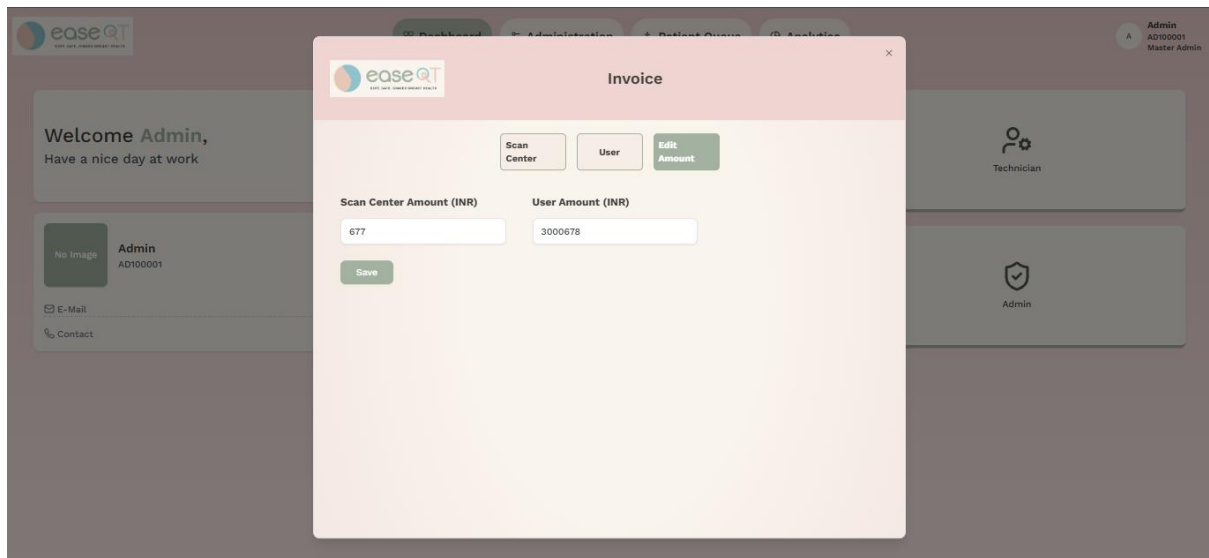


The admin can fix the amount for scan centre and user.


ADMINISTRATION

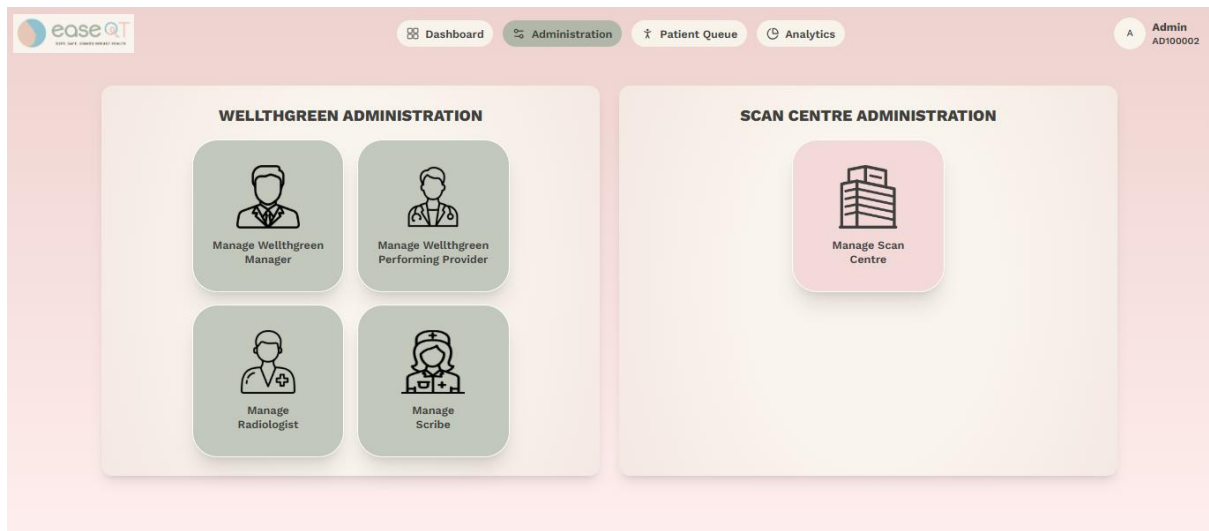

In Administration, the admin can add and manage each user under WellthGreen Administration and scan centres.

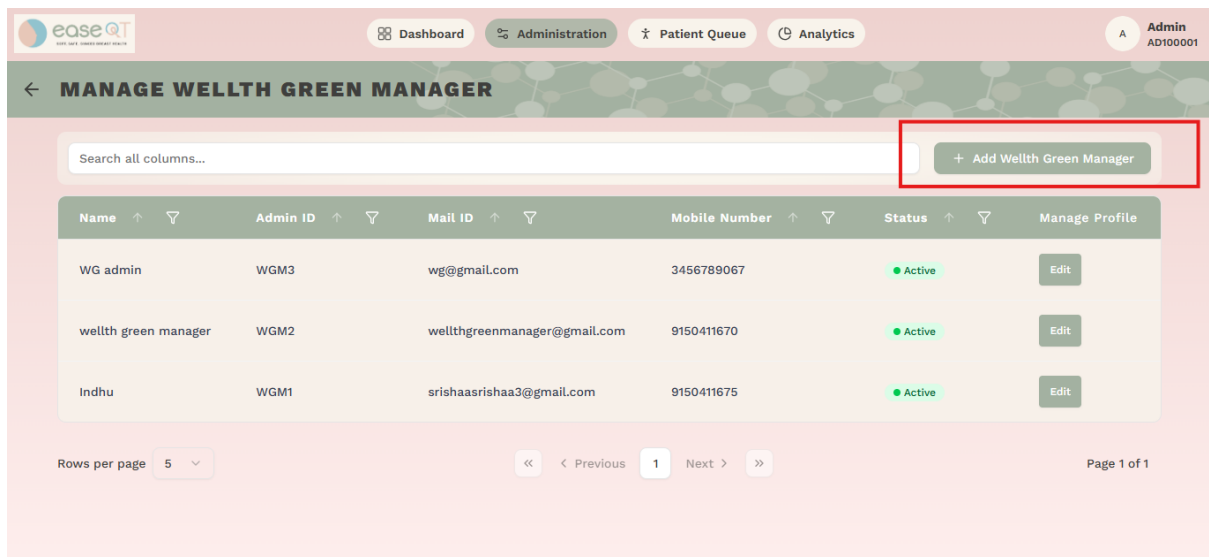


The admin can add Wellthgreen manager by click on 'Add Wellthgreen Manager'.



Dashboard
Administration
Patient Queue
Analytics

A
Admin
AD100001

ADD WELLTH GREEN MANAGER



Add Photo

Personal Details

Full Name *

Contact Number *

IN (+91)

E-Mail *

Date Of Birth *

Aadhar *

PAN *


Driver's License *

Professional Details

Upload Educational Certificates *

Back
Submit

The admin can add the personal and professional details required to add a manager. They can also attach the person's Aadhaar, PAN, driver's license and educational certificates.



Dashboard
Administration
Patient Queue
Analytics

A
Admin
AD100001

MANAGE WELLTH GREEN MANAGER

Name	Admin ID	Mail ID	Mobile Number	Status	Manage Profile
Priyanka	WGM4	priyanka@gmail.com	3456789088	Active	<input type="button" value="Edit"/>
WG admin	WGM3	wg@gmail.com	3456789067	Active	<input type="button" value="Edit"/>
wellth green manager	WGM2	wellthgreenmanager@gmail.com	9150411670	Active	<input type="button" value="Edit"/>
Indhu	WGM1	srishaasrishaa3@gmail.com	9150411675	Active	<input type="button" value="Edit"/>

Rows per page 5

Previous
1
Next

Page 1 of 1

Edit Wellth Green Manager

Manager ID: WGM4

Status: Active

Personal Details

Full Name: Priyanka

Contact Number: 3456789088

E-Mail: priyanka@gmail.com

Date Of Birth: July 1st, 2025

Aadhar Card (.pdf): Choose File No file chosen


PAN Card (.pdf): Choose File No file chosen

Aadhar Document

PAN Document

If the admin wants to change the information in a particular column, they can click the “Edit” option, make the necessary changes, and submit. If they want to deactivate a user, they can change the user’s status from “Active” to “Inactive.”

The same process applies for performing providers, radiologists, and scribes.

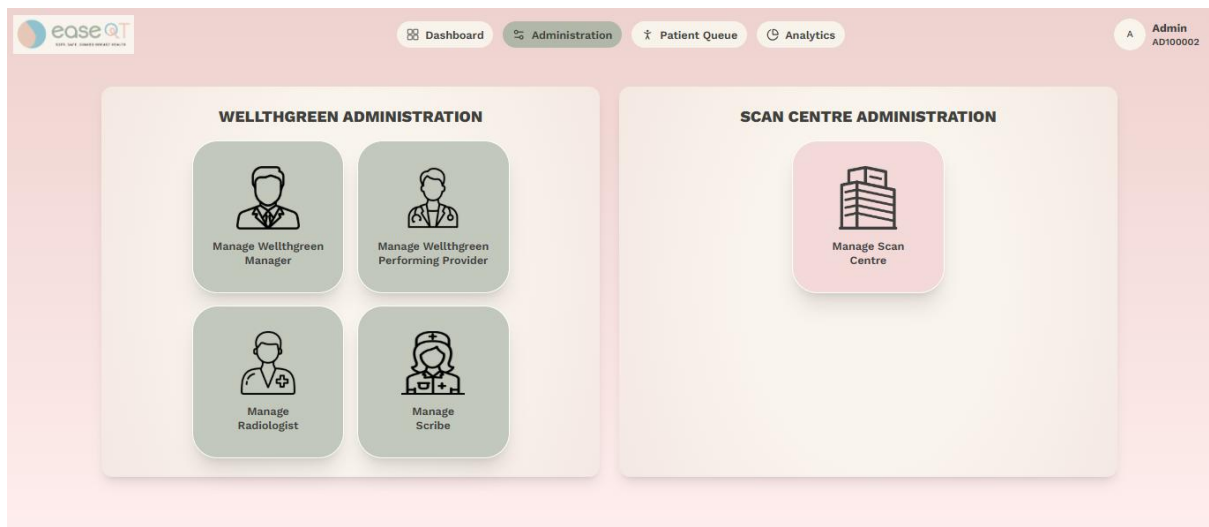

For scan centres, the admin can add new scan centres and assign an admin for each. This admin can then add other users, such as technicians, co-doctors, and performing providers.

PATIENT QUEUE

Dashboard

Administration

Patient Queue

Analytics

Audit Logs

Admin

AD100001

Master Admin

Download Dicom
Download Final Reports
Download Patient Consent
Clear All Filters
Search all columns...

	Date	Scan Centre	Patient ID	Patient Name	Consent	Form	Tech Form	DICOM L	DICOM R	DICOM	Old Report	QT Report	Report Status	Report Delivery	Assign	Pending Remarks	Remarks
<input type="checkbox"/>	Sep 17, 2025	MSC	P100006	Testing12	View	Dc - View	View	(1)	-	View DICOM		View	SIGNED OFF (A)	-	AD10000	I dicom pen...	<div>+ </div>
<input type="checkbox"/>	Sep 18, 2025	MSC	P100006	Testing12	View	S - View	View	(1)	-	View DICOM		View	REVIEWED 2	-	Assign	R1 R2 stuck	<div>+ </div>
<input type="checkbox"/>	Sep 22, 2025	MSC	P100006	Testing12	View	S - View	View	(1)	-	View DICOM		View	SIGNED OFF (A)	-	Assign	No remarks...	<div>+ </div>
<input type="checkbox"/>	Sep 18, 2025	MSC	P100016	sofia	View	S - View	View	(1)	-	View DICOM		View	DRAFT	-	Assign	No remarks...	<div>+ </div>
<input type="checkbox"/>	Sep 19, 2025	MSC	P100016	sofia	View	S - View	View	(1)	-	View DICOM		View	PREDRAFT	-	Assign	No remarks...	<div>+ </div>
<input type="checkbox"/>	Sep 24, 2025	TSC	P100016	sofia	-	Not Filled	Not Yet Started	Not Uploaded	-		-	-	-	-	Assign	No remarks...	<div>+ </div>

Rows per page 50

<<

< Previous


1

Next >

>>

Page 1 of 1

In the patient queue, the admin can view the patient list, patient forms, technician view, and DICOM files. They can fill out reports, check the report status, assign the case to another user, and add pending remarks if needed and add pending remark and can view previous remark also. Click on **View** to proceed with filling out the report. And they can multi select and download DICOM. They can upload and download Patient's old report.



Patient Form

Technician Form

Dicom

Report

General

Right

Left

Impression + Reco

Final Report

← Back

Patient Name: sofia

Patient ID: P100016

Gender: Female

DOB: September 2, 1997 (28 Y)

Patient Form: S Form

Date: 19 Sept 2025

Scan Center: MSC

Other reports

S.No	Report Title	File
No Reports Found		

Other Old reports

S.No	Report Title	File
1	Thermogram	View

Time stamp

Created by	Start	End	Status
sofia	19 Sept	19 Sept	Patient Intake PC

Final Report Preview

Ease QT 10:10 Auto Report

Normal

↺

↻

B

I

U

≡

≡

≡


≡

≡

A

A

📷



NAME	sofia	DOB	September 2, 1997 (28 y)
GENDER	F	SCAN CENTER	MSC, Bangalore
USERID	P100016	DATE OF VISIT	Sep 19, 2025

QT ULTRASOUND BREAST IMAGING

Indication: Screening : Routine Annual Checkup.

History:

Reason for having this QT scan (Patient Stated): ssss.

Current breast symptoms:

- Lump or thickening in the left breast at 6 o'clock since 2 week, size of a pea, right breast at 2 o'clock since 2 week, size of a pea.
- Skin changes in the left breast at 8 o'clock since 2 week, dimpling is noted, right breast at 3 o'clock since 2 week, dimpling is noted.

In the report section, the admin can review and edit reports completed by the radiologist. There are “Reviewed 1 Correct” and “Reviewed 2 Edit” buttons available for the admin. If the radiologist has completed the report correctly, the admin will select “Reviewed 1 Correct.” If it is incorrect, the admin will select “Reviewed 2 Edit.” After verifying the report, the admin will sign off on it.

[← Back](#)
 Patient Name: sofia
 Patient ID: P100016
 Gender: Female
 DOB: September 2, 1997 (28 Y)
 Patient Form: S Form
 Date: 19 Sept 2025
 Scan Center: MSC

☒ Right Breast

Impression		Recommendation	
0	Incomplete/Indeterminate Study	0	Repeat QT Imaging
1	Likely Benign Findings	1	Routine Annual Screening
1a	Likely Benign Findings + Architectural Distortion	1a	Routine Annual Screening
2	Low-Risk Characteristics	2	Short-Term Follow-Up 6m
2a	Low risk - too small	2a	Short-term Follow-up 6m
3	Indeterminate Findings: Non palpable mass	3	Growth Rate Assessment 3m
3a	Indeterminate: Not Definitely Abnormal	3a	Confirm Benign Nature with Targ USG
3b	Indeterminate: Calcifications	3b	Targ USG for Cystic Mass Evaluation
3c	Indeterminate: Cysts	3c	Clinical consultation with a breast specialist or health care provider for Calcification Management
3d	Indeterminate : Postop changes + fluid/ seroma	3d	Clinical Consultation
3e	Indeterminate: Post-surgical vs. Recurrence	3e	Clinical consultation + Surveillance
3f	Indeterminate: Likely Fibrosis	3f	Clinical Consultation + Surveillance
3g	Indeterminate: Long Standing Findings	3g	If persists - Targ USg
4	Pain But No LUMP + QT Finding	4	Clinical Monitoring
4a	Palp Lump - No QT Finding or Ab/Normal image findings	4a	Breast specialist Evaluation, neg never preclude/Targ USG.
4b	Palpable Lump + QT Finding		

Other reports

S.No	Report Title	File
No Reports Found		

Other Old reports

S.No	Report Title	File
1	Thermogram	View

Time stamp

Created by	Start	End	Status
sofia	19 Sept 2025	19 Sept 2025	Patient Intake Fc

The admin can give the impression and recommendation by clicking on it. They can add extra impression and recommendation.

Time stamp

Created by	Start	End	Status
sofia	19 Sept 2025 18:29	19 Sept 2025 06:09	Patient Intake Form Filled
Technician	19 Sept 2025 06:10	19 Sept 2025 06:14	Technology Form Filled
Wellthgreen	19 Sept 2025 06:15	19 Sept 2025 06:16	Predraft

Status saver

Load Template

Reset to Default

Predraft

Draft

Reviewed 1 Correct

Reviewed 1 Edit

Reviewed 2 Correct

Reviewed 2 Edit

Insert Signature

Sign Off

Final Report Preview

Ease QT 10:10 Auto Report

Normal

↶

↷

B

I

U

☰

☷

☱

☲

≡

A

🗑️

🖨️

breast at 2 o'clock since 2 week, size of a pea.

- Skin changes in the left breast at 8 o'clock since 2 week, dimpling is noted, right breast at 3 o'clock since 2 week, dimpling is noted.

TECHNIQUE: Transmission and reflection multiplanar 3-dimensional ultrasound imaging of both breasts was performed using the QT Ultrasound Series 2000 Model-A scanner. Breast density was determined using the Quantitative Breast Density calculator. Images were reviewed in the QTviewer v2.6.2. The nipple-areolar complex, skin, Cooper's ligaments, breast fat distribution, penetrating arteries and veins, glandular and ductal tissues were evaluated. Images were reviewed in coronal, transaxial and sagittal planes.

BREAST IMPLANTS:

The QT scan shows implant which are subpectoral in position, with speed of sound consistent with . The displacement is noted in right side.

SYMMETRY:

The breast examination indicates symmetrical size and shape.

RIGHT BREAST FINDINGS:

Density: The right breast tissue is heterogeneously dense with a fibroglandular ratio of 4:6.

There is an option called "Reset to Default." If the report was done incorrectly, they can use this option to revert it to the default state.

Patient Form

Technician Form

Dicom

Report

General

Right

Left

Impression + Reco

Final Report

← Back

Patient Name: sofia
Patient ID: P100016
Gender: Female
DOB: September 2, 1997 (28 Y)
Patient Form: S Form
Date: 19 Sept 2025
Scan Center: MSC

Other reports

S.No	Report Title	File
No Reports Found		

Other Old reports

S.No	Report Title	File
1	Thermogram	View

Time stamp

Created by	Start	End	Status
sofia	19 Sept 2025 18:29	19 Sept 2025 06:09	Patient Intake Form Filled Completed
Technician	19 Sept 2025 06:10	19 Sept 2025 06:14	Technologist Form Filled
Wellthgreen Scribe	19 Sept 2025 06:14	19 Sept 2025	Prodraft

Final Report Preview

Ease QT 10.10 Auto Report

Normal

B I U

NAME

sofia

DOB

September 2, 1997 (28 Y)

GENDER

F

SCAN CENTER

MSC, Bangalore

USERID

P100016

DATE OF VISIT

Sep 19, 2025

QT ULTRASOUND BREAST IMAGING

Indication: Screening : Routine Annual Checkup.

History:

Reason for having this QT scan (Patient Stated): ssss.

Current breast symptoms:

- Lump or thickening in the left breast at 9 o'clock since 2 week, size of a pea, right breast at 2 o'clock since 2 week, size of a pea
- Skin changes in the left breast at 8 o'clock since 2 week, dimpling is noted, right breast at 3 o'clock since 2 week, dimpling is noted.

TECHNIQUE: Transmission and reflection multipolar 3-dimensional ultrasound imaging of both breasts was performed using the QT Ultrasound Series 2000 Model-A-scanner. Breast density was determined using the Quantitative Breast Density calculator. Images were reviewed in the QTviewer v2.0.2. The nipple-areolar complex, skin, Cooper's ligaments, breast fat distribution, penetrating arteries and veins, glandular and ductal tissues were evaluated. Images were reviewed in coronal, transaxial and sagittal planes.

The admin can view the full report here. If they want to make changes or load a template, **they need to ensure the easeqt 10.10 is off**. If they want to access the report, they need to switch on the 10.10 auto-report feature.

AUDIT LOGS

Dashboard

Administration

Patient Queue

Analytics

Audit Logs

Admin

AD100001

Master Admin

User ID	Data	Timestamp
AD100001	Logged In Successfully	20/9/2025, 3:39:33 pm
WGD5	Logged In Successfully	20/9/2025, 3:32:24 pm
CMSCR1	Logged In Successfully	20/9/2025, 3:20:17 pm
U	Logged In Successfully	20/9/2025, 2:55:06 pm
U	Logged In Successfully	20/9/2025, 2:31:44 pm
AD100003	Logged In Successfully	20/9/2025, 2:28:45 pm
CMSCR1	Logged In Successfully	20/9/2025, 2:27:32 pm
AD100001	Logged In Successfully	20/9/2025, 2:25:52 pm
AD100003	Logged In Successfully	20/9/2025, 1:34:18 pm
AD100001	Logged In Successfully	20/9/2025, 1:33:46 pm
AD100003	Logged In Successfully	20/9/2025, 5:40:02 am

Rows per page

50

<<

< Previous

1

Next >

>>

Page 1 of 20



Copyright © Wellthgreen

The admin can view all logged in activity of users and patients in audit log.

ANALYTICS



The admin can track analytics for all scan centres, users and their own activities. For their personal analytics, they can view the total count, number of edits, number of correct reports, counts for the previous six months, TAT (Turnaround Time), scan indications, artifacts and recommendation codes. The admin can download as an excel of all users and scan centre.

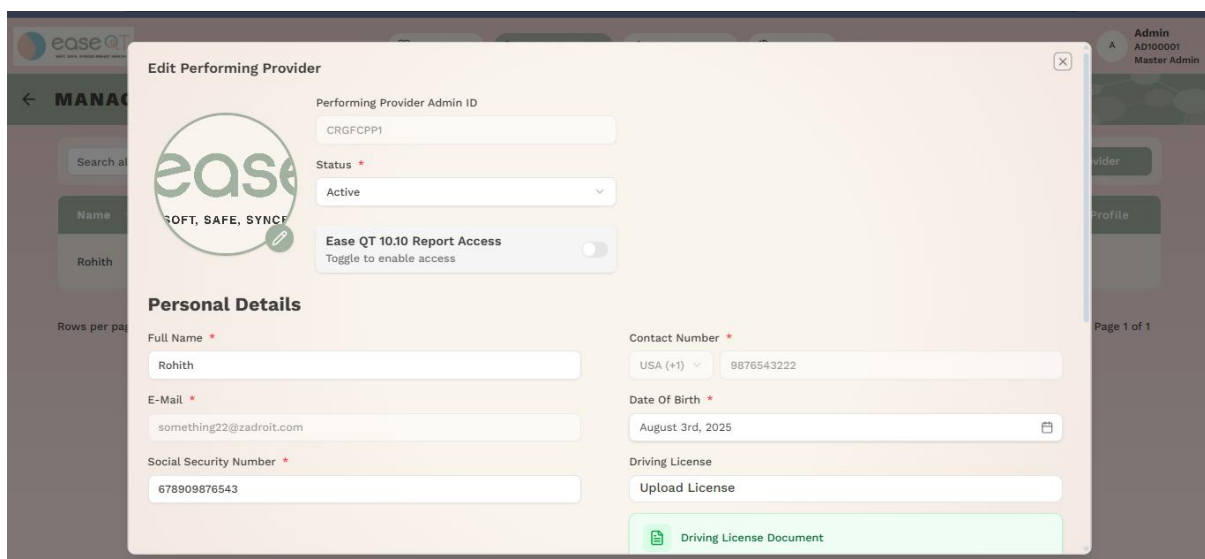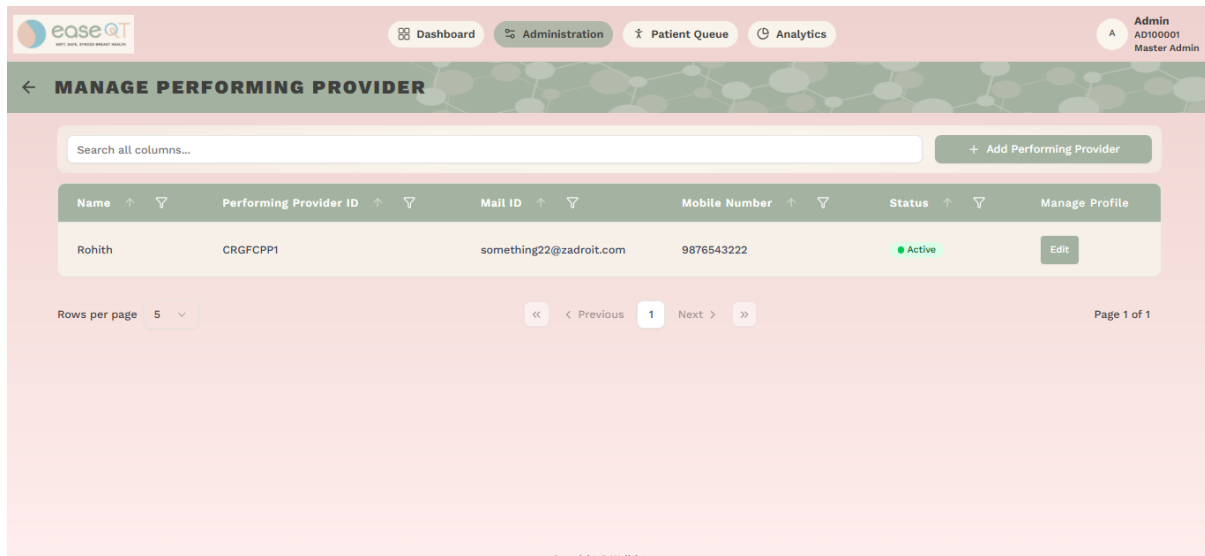

If the performing provider and reviewer want to view the general details, left breast details, and right breast details, the admin need to click edit their profile and enable the Ease QT 10.10 report access.

GENERAL TIPS

1. If you open an item in edit mode, make sure to save it before exiting; otherwise, others won't be able to access it.
2. After addressing pending remarks, add a new remark saying 'Done' and set 'Assign to' as None."
3. The Filter option is available for sorting.
4. For bulk downloads, use the checkboxes to select patients and then download in bulk.
5. To add a pending remark, click the + icon.
6. When viewing a report, always click Back to exit-otherwise, the case you opened may freeze.
7. Use Image Preview to add images to a report.
8. To refresh the page, press Ctrl + Shift + R.
9. For DICOM downloads and voice dictation, use Chrome or Edge. (For Firefox, click Save. Do not use Internet Explorer.)
10. Always log out properly-do not just close the browser or tab.